



Thermo Scientific™ General Purpose Centrifuge Series

To select the right centrifuge for your lab,
choose an option at the top of the screen,
or click below to see what's new

Let's get started

ДИА•М
современная лаборатория

www.dia-m.ru
заказ on-line





Continuous innovation

New industrial design with improved ergonomics for ease of use and touch screen capabilities to help save you time

**View
GUI**



State-of-the-art touchscreen interface

**View
safe operation**



Safety and ergonomic enhancements

**View
enhancements**



Enhancements to make your day in the lab easier



Graphic User Interface (GUI)

Modern, intuitive touchscreen interface




User-intuitive screen enables easy access and programming for faster results

[View](#)



New Lighthouse Mode provides run status visibility from across the lab, now with adjustable-volume alerts

[View](#)

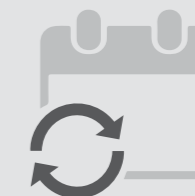





Scheduling features and preventative maintenance queues to help reduce downtime


[View](#)



Simple touchscreen programming

Click below to discover the new touchscreen features

- Easily accessible **menu bar** 
- Keep an eye on **centrifuge health** 
- Time-saving **pre-temp button** 
- Customize **quick spins** 
- Access **helpful info** with one touch 

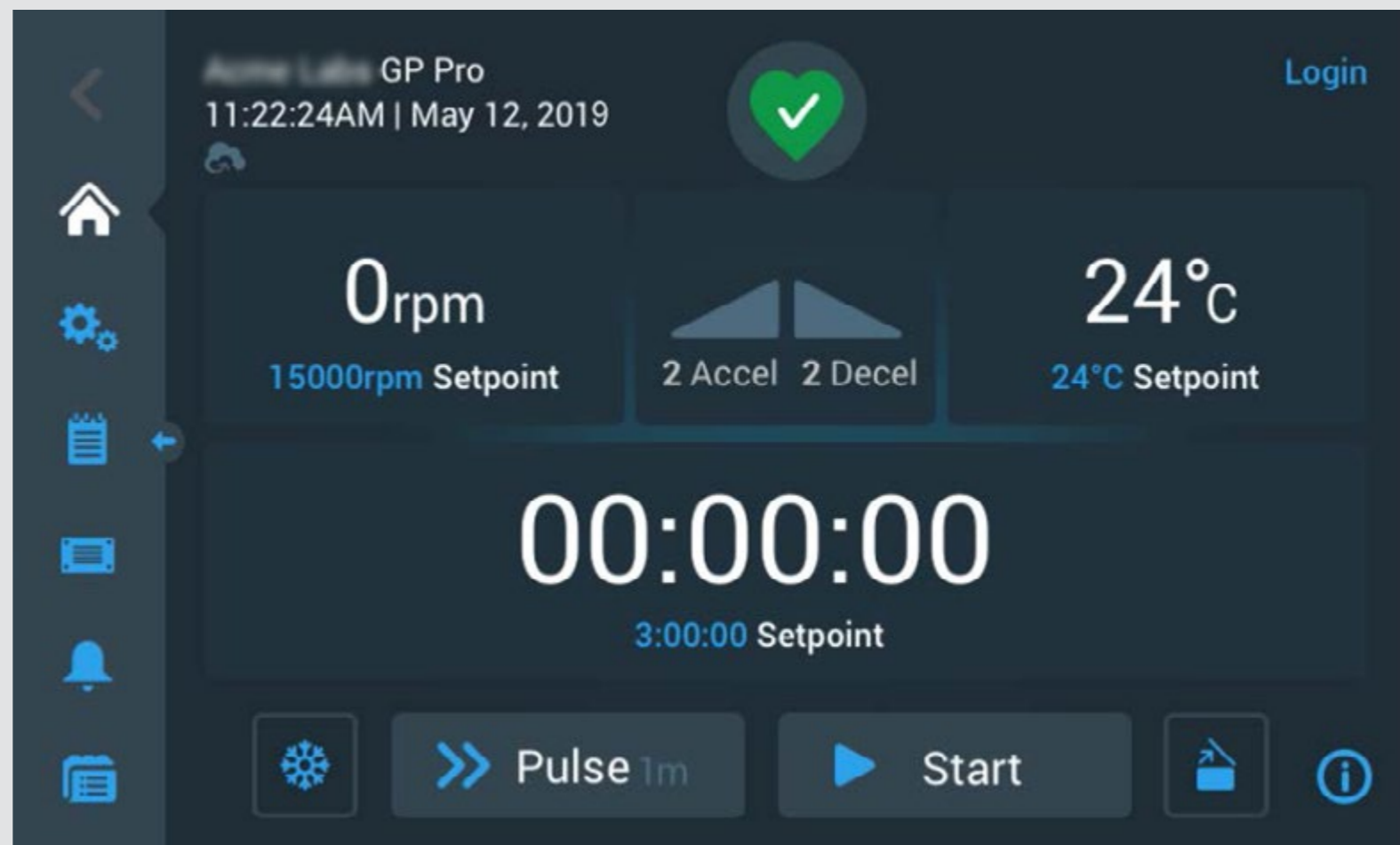




Lighthouse mode



Track progress without having to step away from your work

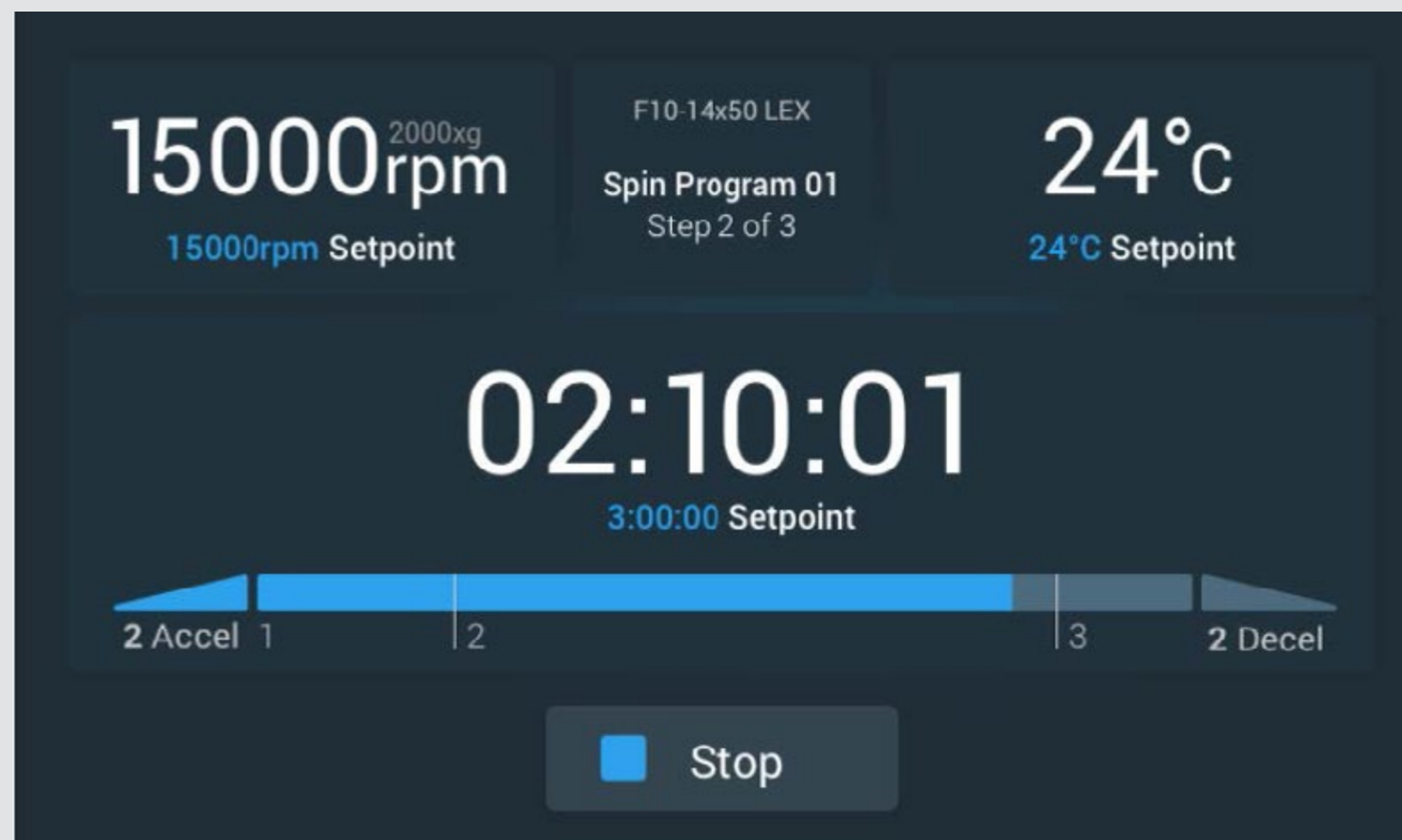
- Easily **monitor your run progress** from across the lab with enlarged time and tracker bar
- Plus, **end-of-run alerts** can be programmed with customized volume settings





Smart controls

- Pre-programmed start-up of compressor for temperature-ready centrifuge as soon as you arrive in the lab 
- Shut-down programming can extend compressor life and reduce energy costs 

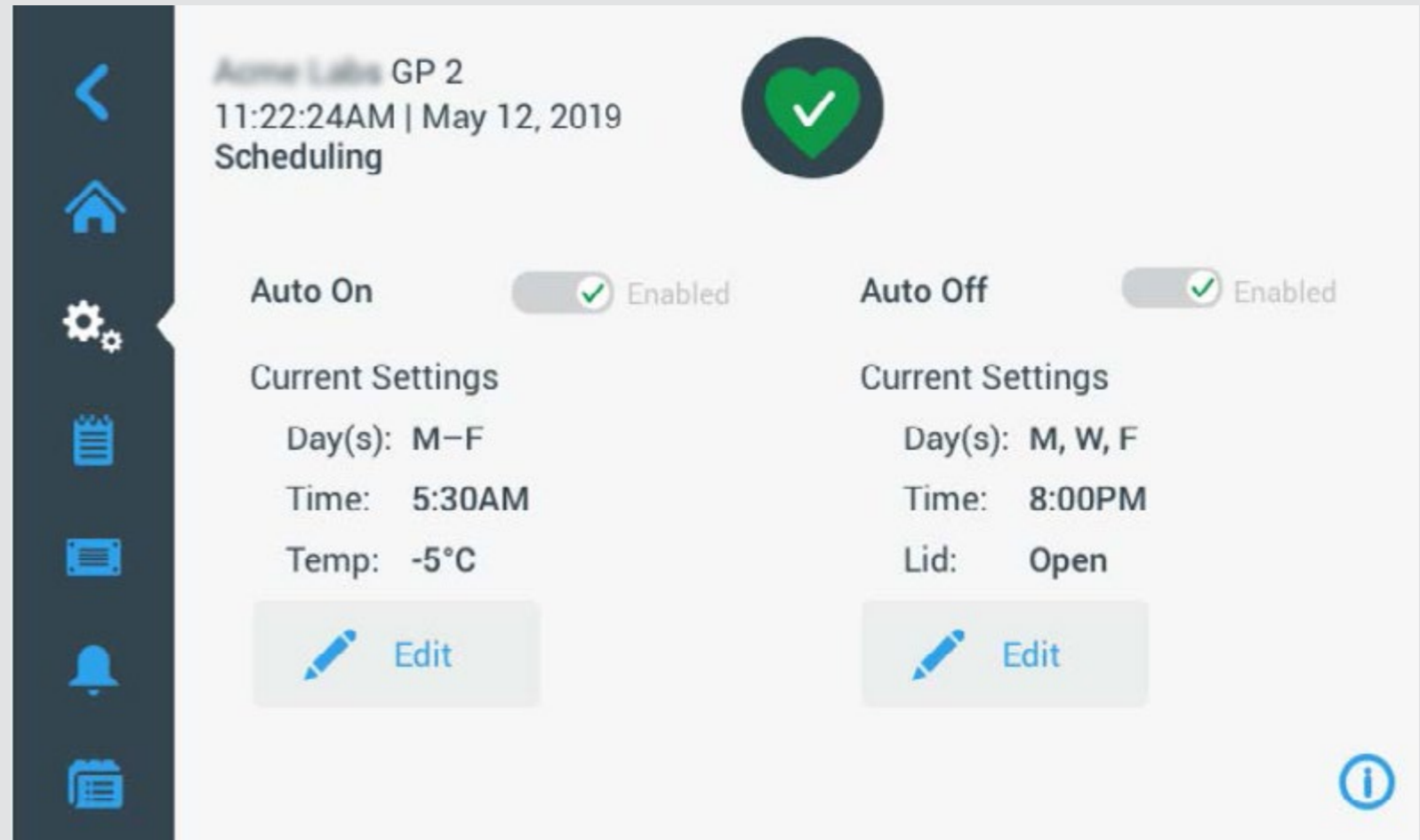

Acme Labs GP 2
11:22:24AM | May 12, 2019
Scheduling

Auto On Enabled

Auto Off Enabled

Current Settings


Day(s): M-F
Time: 5:30AM
Temp: -5°C

Edit

Current Settings

Day(s): M, W, F
Time: 8:00PM
Lid: Open

Edit



Safety features

Improved reach-in during rotor installation for safer loading and unloading

Hand grip cue on the lid indicates the easiest and best one-handed closure location

Reduced pinch-points designed to prevent accidental injury. Motorized lid latch with safety protection

[View](#)

[View](#)





Enhancements

Accessory tray holds tools or samples as you prepare for your run

[View](#)

Front USB port for easier plug in and data transfer

[View](#)

Other new benefits include:

- compliance with most-recent applicable regulatory and safety standards
- available energy-saving settings
- over 20 rotor options to choose from



Consistent performance

Thermo Scientific technologies designed for solid performance and consistent results in a wide range of applications

View performance



Broad choice of speed, capacity and application

View Auto-Lock™



Simple, fast rotor exchange

View Fiberlite™



Outstanding rotor ergonomics and performance

View ClickSeal™




One-handed biocontainment closures



Find efficiency with more than 20 rotor options

Achieve faster results with the ability to hold up to 196 blood tubes, 40 x 50 mL conical tubes, or 24 microplates in one run

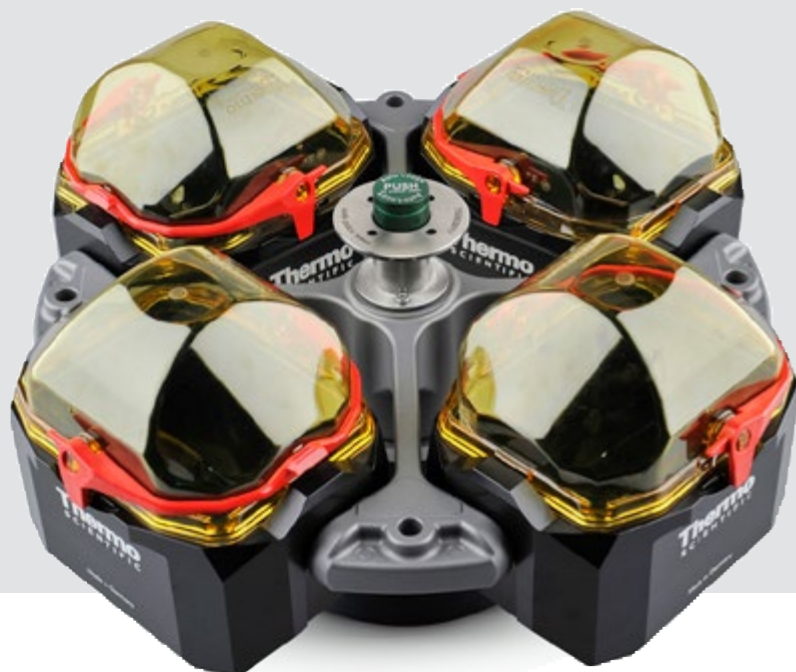

- With fixed-angle **rotor speeds over 25000 x g**, and swinging buckets rotor speeds over 6000 x g, even the toughest separations can be completed better and faster >
- **Motorized lid latch** reduces closure time and is easy to operate/handle >
- **Slide-coat finish** on swinging buckets rotors reduces the need for lubrication, for more consistent separations and less chance of grease contamination >



Auto-Lock™ Rotor Exchange

Helps save time. Faster answers.

- Secure, ergonomic push-button rotor exchange in as little as 3 seconds
- Trouble-free installation and removal with **no tools required**
- **Full chamber access** for cleaning convenience and minimizes risk of drip tray contamination
- **Flexibility to quickly switch between applications** and evolve with the changing needs of your laboratory



Click below to learn how Fiberlite™ Rotors let researchers focus on their work – worry-free

1

Biocontainment



2

Design



3

Performance



4

Ergonomics



5

Value



Fiberlite™ Carbon Fiber Rotors

Outstanding ergonomics and performance

Lighter materials, advanced designs and industry experience all contribute to the development of superior rotor products with better performance, efficiency and built-in safety features



5 Ways

Fiberlite™ Rotors let researchers focus on their work – worry-free

1**2****3****4****5**

1 Biocontainment

- Fiberlite™ Rotor / lid combinations, and liquid containment annulus offer features that protect the user, the centrifuge and the samples from hazardous contamination
- Biocontainment certified by Public Health England, Porton Down, UK



5 Ways

Fiberlite™ Rotors let researchers focus on their work – worry-free

1

2


3

4

5



Innovative molded-design to protect tubes from cracking, sample loss, and centrifuge damage.

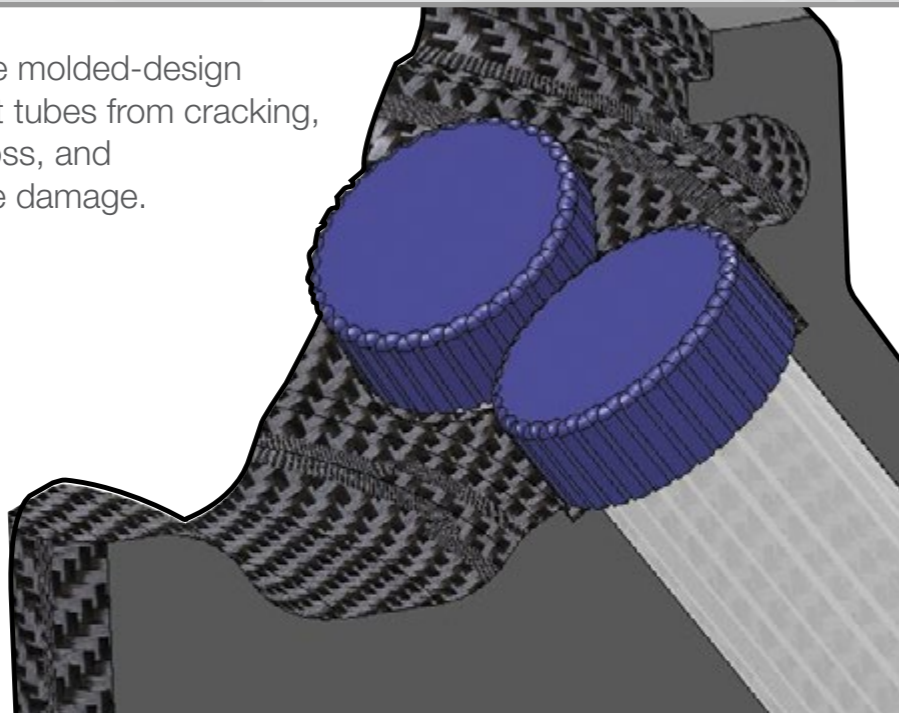

2 Design

- Corrosion-resistance and fatigue-resistance maintain the rotor's structural integrity for outstanding durability
- Fiberlite™ Rotors are designed to “unwind” rather than split or break, so your centrifuge is protected from damage in the unlikely event of a rotor failure
- Carbon fiber offers a unique ability to mold each cavity to perfectly fit tubes, bottles and other consumables, enabling faster spins for tougher separations



5 Ways

Fiberlite™ Rotors let researchers focus on their work – worry-free

1

2

3

4

5



Time savings with carbon fiber rotors

Acceleration and deceleration rate comparison of 6 x 250 mL capacity floor model rotors*

Fiberlite™ Carbon Fiber Rotor	1:35/1:15	Savings (2:55 minutes)
Aluminum rotor	3:45/2:00	

Accel/decel rates

3 Performance

- Thermo Scientific™ Fiberlite™ Rotors are each balanced and tested before they leave the factory for precision performance
- Automatic rotor recognition eliminates the need for manual input**
- Fiberlite™ Rotors can reach higher speeds than metal equivalents, improving separations in tubes such as 15 or 50 mL conicals

* Based on a comparison with manufacturers' published specifications.

**After reaching speeds of 200 rpm; please see manual for details.



5 Ways

Fiberlite™ Rotors let researchers focus on their work – worry-free

1

2

3


4

5



Weight savings with carbon fiber rotors

Weight comparison of fully loaded 6 x 250 mL capacity floor model rotors*

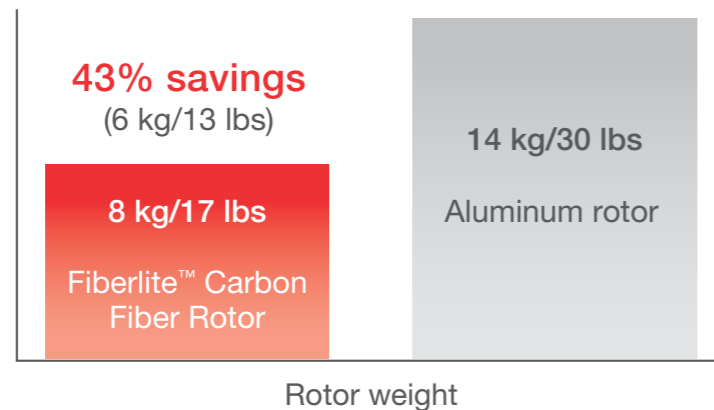

4 Ergonomics

- Auto-Lock™ allows push-button rotor exchange in as little as 3 seconds with no tools or complicated procedures required
- Carbon fiber is lighter in weight allowing for easier lifting in and out of the centrifuge
- Unique integrated handle design allows for easy gripping to exchange rotors for a variety of applications
- Easy rotor exchange simplifies cleaning and disinfecting procedures

* Based on a comparison with manufacturers' published specifications.



5 Ways


Fiberlite™ Rotors let researchers focus on their work – worry-free

- 1
- 2
- 3
- 4
- 5



Warranty with carbon fiber rotors

Average warranty periods for metal rotors compared with Fiberlite™ Carbon Fiber Rotors*

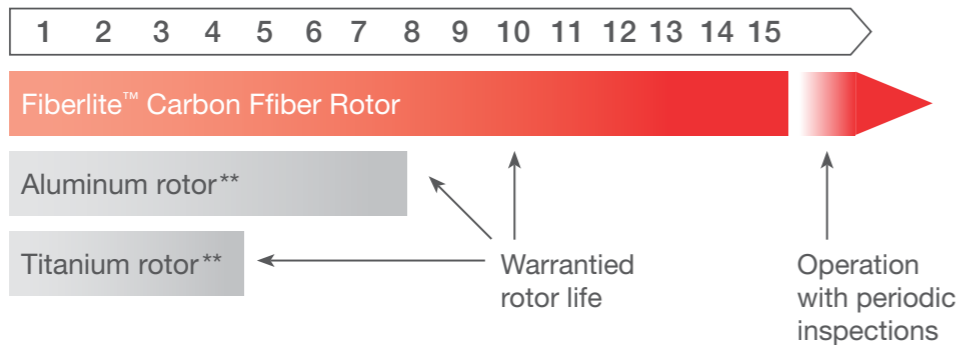

5 Value

- On-site clinics and maintenance are available to ensure proper use and care, extending the life of your rotors
- Our continuously improving rotor designs have increased the capacities available on the same size centrifuge, resulting in outstanding efficiency and throughput across the range
- Warranty to support fewer replacements

* Average warranty periods were calculated based on industry average of years an aluminum or titanium rotor may be covered under warranty per manufacturers' published specifications.

**Warranty coverage may vary by rotor. Please refer to manufacturer for specific warranty coverage for each rotor.



ClickSeal™ Biocontainment Lids

One-handed sample protection

- Biocontainment sealing options for glove-friendly, one-handed open/close capability are ergonomically designed for both right- and left-handed operations
- Simple operation for all laboratory users, eliminating multi-turn screw caps and complicated high pressure clips
- Biocontainment certification by Public Health England, Porton Down, UK





Meet our family of centrifuges

Click on the sizing tool below to explore the best options for your lab space and applications

General Purpose

1.6 L Ventilated [➤](#)

General Purpose

1.6 L Refrigerated [➤](#)

General Purpose

4 L Floor Model

Ventilated &

Refrigerated [➤](#)

General Purpose

4 L Ventilated [➤](#)

General Purpose

4 L Refrigerated [➤](#)

View All [➤](#)



Specifications and ordering information

Specifications	Multifuge X4 Pro		Multifuge X4R Pro		Megafuge ST4 Plus		Megafuge ST4R Plus	
	IVD/MD	General Use Only	IVD/MD	General Use Only	IVD/MD	General Use Only	IVD/MD	General Use Only
Swinging Buckets Rotors								
Maximum Capacity	4 x 1000 mL with TX-1000 rotor							
Maximum Speed	6000 rpm with BIOShield 1000A rotor				5300 rpm with BIOShield 1000A rotor			
Maximum RCF	7164 x g with BIOShield 1000A rotor				5590 x g with BIOShield 1000A rotor			
Fixed Angle Rotors								
Maximum Capacity	6 x 250 mL with F14-6x250LE rotor				6 x 100 mL with F15-6x100 or F10-6x100 rotor			
Maximum Speed	15200 rpm with F21-48x2, Microliter 48x2 or Microliter 30x2							
Maximum RCF	25830 x g with Microliter 30x2							
Centrifuge Details								
Control System	Microprocessor							
Drive System	Direct, Brushless induction low profile motor							
Accel/Decel Profiles	Choices of 9 Accel / 10 Decel Rates							
Material	front, top, sides: PC-ABS plastic painted with 2K PUR-based lacquer; metal housing: cold rolled steel with powder paint; centrifuge bowl: stainless steel; centrifuge feet: rubber							
Safety Features	SMARTSpin imbalance detection, finger-pinch prevention, crash-proof construction							
Warranty	5 Year Limited							
Includes	printed English user manual, USB stick with 28 languages of user manual, corrosion resistance oil (70009824), power cord, declaration of conformity							
Noise Level	<64 dB for refrigerated and <68 dB for ventilated							
Temperature Range	N/A		-10°C to +40°C		N/A		-10°C to +40°C	
Number of Compatible Rotors	17		17		12		12	
Display Details								
Display	Full color touch screen interface				Push button and high contrast LCD interface			
Run Time	max 99h, 59min, 59sec; also with continuous and customizable pulse				max 9h, 59min; also with continuous and pulse modes			
Timer	count up or down				count down only			
Program Storage	up to 100 run programs				up to 6 run programs			
Programmability	via easy access touch screen menu				via 6 direct access buttons			
Dimensions, Weight								
Width (inches) Exterior	22.3		29.4		22.3		29.4	
Width (cm) Exterior	56.6		74.6		56.6		74.6	
Depth (inches) Exterior	27.2		27.2		27.2		27.2	
Depth (cm) Exterior	69.0		69.0		69.0		69.0	
Height (inches) Exterior	14.3		14.2		14.3		14.2	
Height (cm) Exterior	36.2		36.1		36.2		36.1	
Height (inches) Open Cover	33.9		33.9		33.9		33.9	
Height (cm) Open Cover	86.1		86.0		86.1		86.0	
Footprint Dimensions (inches)	22.3 x 23.1		29.4 x 23.0		22.3 x 23.1		29.4 x 23.0	
Footprint Dimensions (cm)	56.6 x 58.8		74.6 x 58.4		56.6 x 58.8		74.6 x 58.4	
Net Weight (lbs)	196		265		196		265	
Net Weight (kgs)	89		120		89		120	
Certifications and Standards								
IVD/MD	X		X		X		X	
Certifications	ENUL/IEC 61010-1, ENUL/IEC 61010-2-020 and 61010-2-101, ENUL/IEC 61326-1, ENUL/IEC 61326-2-6, FCC part 15							
Standards	ISO 14971, ISO 13485, ISO 9001							
Ordering Information								
	Multifuge X4 Pro		Multifuge X4R Pro		Megafuge ST4 Plus		Megafuge ST4R Plus	
	IVD/MD	General Use Only	IVD/MD	General Use Only	IVD/MD	General Use Only	IVD/MD	General Use Only
230V [§]	75009500*	75009900	75009515* / 75009615*	75009915* / 75009815*	75009503*	75009903	75009518* / 75009618*	75009918* / 75009818*
120V	75009501*	N/A	75009516*	N/A	75009504*	N/A	75009519*	N/A

§ 220V / 60 Hz [Korea, etc.]

*Indicates a centrifuge, rotor and/or adapter product approved for use in clinical applications; see user manual for details.



One part number. The complete package.

Order your centrifuge, rotor and accessories with a single part number designed for your application

View

Packages with swinging buckets rotors developed for 50 and 15 mL conical tubes. Excellent for cell culture and general life science applications with up to 96 x 15 mL conical tube capacities

View

With up to 196 standard 5/7 mL blood tube capacities in a benchtop footprint, process more samples for faster patient or study results





Cell culture packages

Centrifuge	Package PN	Package Includes
4L High Performance including Touch Screen User Interface		
40 x 50 mL conical tubes up to 4200 rpm/4122 x g (R) or 3800 rpm/3374 x g (V)		
Multifuge X4R Pro-MD, 120V	75016045*	TX-1000 rotor with Auto-Lock™ Buckets (set/4) ClickSeal™ biocontainment lids (set/4) 50 mL conical tube adapters (set/4, 40 tube capacity) 15 mL conical tube adapters (set/4, 96 tube capacity)
Multifuge X4R Pro, 208-240V	75016066	
Multifuge X4 Pro-MD, 120V	75016046*	
Multifuge X4 Pro, 208-240V	75016067	

Centrifuge	Package PN	Package Includes
4L High Performance including Touch Screen User Interface		
28 x 50 mL conical tubes up to 4700 rpm/4816 x g		
Multifuge X4 Pro, 208-240V	75016077	TX-750 rotor with Auto-Lock Round Buckets (set/4) ClickSeal biocontainment lids (set/4) 50 mL conical tube adapters (set/4, 28 tube capacity)

Centrifuge	Package PN	Package Includes
4L Standard Performance		
40 x 50 mL conical tubes up to 4200 rpm/4122 x g (R) or 3800 rpm/3374 x g (V)		
Megafuge ST4R Plus-MD, 120V	75016049*	TX-1000 rotor with Auto-Lock Buckets (set/4) ClickSeal biocontainment lids (set/4) 50 mL conical tube adapters (set/4, 40 tube capacity) 15 mL conical tube adapters (set/4, 96 tube capacity)
Megafuge ST4R Plus, 208-240V	75016069	
Megafuge ST4 Plus-MD, 120V	75016050*	
Megafuge ST4 Plus, 208-240V	75016070	

Centrifuge	Package PN	Package Includes
4L Standard Performance		
28 x 50 mL conical tubes up to 4700 rpm/4816 x g		
Megafuge ST4R Plus, 208-240V	75016083	TX-750 rotor with Auto-Lock Round Buckets (set/4) ClickSeal biocontainment lids (set/4) 50 mL conical tube adapters (set/4, 28 tube capacity) 15 mL conical tube adapters (set/4, 56 tube capacity)

*Indicates a centrifuge, rotor and/or adapter product approved for use in clinical applications; see user manual for details.



Rotor-only cell culture packages

Rotor	Package PN	Package Includes
TX-1000 Conical Tube Package		
TX-1000 Swinging Buckets Rotor	75016105*	TX-1000 rotor cross with Auto-Lock™ Buckets (set/4) ClickSeal™ biocontainment lids (set/4) 50 mL conical tube adapters (set/4, 40 tube capacity) 15 mL conical tube adapters (set/4, 96 tube capacity)
TX-1000 Microplate Package		
TX-1000 Swinging Buckets Rotor	75016104*	TX-1000 rotor cross with Auto-Lock Buckets (set/4) ClickSeal biocontainment lids (set/4) Microplate adapters (set/4; 24 plate capacity, 8 deepwell capacity)
TX-750 Conical Tube Package		
TX-750 Swinging Buckets Rotor	75016110*	TX-750 rotor cross with Auto-Lock Round Buckets (set/4) ClickSeal biocontainment lids (set/4) 50 mL conical tube adapters (set/4, 28 tube capacity) 15 mL conical tube adapters (set/4, 56 tube capacity)
TX-400 Conical Tube Package		
TX-400 Swinging Buckets Rotor	75016108*	TX-400 rotor cross with Auto-Lock Buckets (set/4) ClickSeal biocontainment lids (set/4) 50 mL conical tube adapters (set/4, 16 tube capacity) 15 mL conical tube adapters (set/4, 36 tube capacity)
TX-200 Conical Tube Package		
TX-200 Swinging Buckets Rotor	75016106*	TX-200 rotor cross with Auto-Lock Buckets (set/4) ClickSeal biocontainment lids (set/4) 50 mL conical tube adapters (set/4, 4 tube capacity) 15 mL conical tube adapters (set/4, 20 tube capacity)

*Indicates a centrifuge, rotor and/or adapter product approved for use in clinical applications; see user manual for details.



Clinical packages

Centrifuge	Package PN	Package Includes
4L High Performance including Touch Screen User Interface		
196 x 5/7 mL tubes up to 4200 rpm/4122 x g (R) or 3800 rpm/3374 x g (V)		
Multifuge X4R Pro-MD, 120V	75016033*	TX-1000 rotor with Auto-Lock™ Buckets (set/4)
Multifuge X4R Pro-MD, 208-240V	75016057*	ClickSeal™ biocontainment lids (set/4)
Multifuge X4 Pro-MD, 120V	75016034*	5/7 mL blood tube adapters (set/4, 196 tube capacity)
Multifuge X4 Pro-MD, 208-240V	75016058*	10 mL blood tube adapters (set/4, 148 tube capacity)

Centrifuge	Package PN	Package Includes
4L High Performance including Touch Screen User Interface		
148 x 5/7 mL tubes up to 4700 rpm/4816 x g		
Multifuge X4R Pro-MD, 208-240V	75016081*	TX-750 rotor with Auto-Lock Round Buckets (set/4) ClickSeal biocontainment lids (set/4) 5/7 mL blood tube adapters (set/4, 148 tube capacity) 10 mL blood tube adapters (set/4, 84 tube capacity)



Centrifuge	Package PN	Package Includes
4L Standard Performance		
196 x 5/7 mL tubes up to 4200 rpm/4122 x g (R) or 3800 rpm/3374 x g (V)		
Megafuge ST4R Plus-MD, 120V	75016037*	TX-1000 rotor with Auto-Lock Buckets (set/4)
Megafuge ST4R Plus-MD, 208-240V	75016060*	ClickSeal biocontainment lids (set/4)
Megafuge ST4 Plus-MD, 120V	75016038*	5/7 mL blood tube adapters (set/4, 196 tube capacity) 10 mL blood tube adapters (set/4, 148 tube capacity)

Centrifuge	Package PN	Package Includes
4L Standard Performance		
148 x 5/7 mL tubes up to 4700 rpm/4816 x g		
Megafuge ST4R Plus-MD, 208-240V	75016079*	TX-750 rotor with Auto-Lock Round Buckets (set/4) ClickSeal biocontainment lids (set/4) 5/7 mL blood tube adapters (set/4, 148 tube capacity) 10 mL blood tube adapters (set/4, 84 tube capacity)

*Indicates a centrifuge, rotor and/or adapter product approved for use in clinical applications; see user manual for details.



Rotor-only clinical packages

Rotor	Package PN	Package Includes
TX-1000 Blood Tube Package		
TX-1000 Swinging Buckets Rotor	75016103*	TX-1000 rotor cross with Auto-Lock™ Buckets (set/4) ClickSeal™ biocontainment lids (set/4) 5/7 mL blood tube adapters (set/4, 148 tube capacity) 10 mL blood tube adapters (set/4, 84 tube capacity)

Rotor	Package PN	Package Includes
TX-750 Blood Tube Package		
TX-750 Swinging Buckets Rotor	75016109*	TX-750 rotor cross with Auto-Lock Round Buckets (set/4) ClickSeal biocontainment lids (set/4) 5/7 mL blood tube adapters (set/4, 148 tube capacity) 10 mL blood tube adapters (set/4, 84 tube capacity)

Rotor	Package PN	Package Includes
TX-400 Blood Tube Package		
TX-400 Swinging Buckets Rotor	75016107*	TX-400 rotor cross with Auto-Lock Buckets (set/4) ClickSeal biocontainment lids (set/4) 5/7 mL blood tube adapters (set/4, 76 tube capacity) 10 mL blood tube adapters (set/4, 56 tube capacity)

Rotor	Package PN	Package Includes
TX-200 Blood Tube Package		
TX-200 Swinging Buckets Rotor	75016111*	TX-200 rotor cross with Auto-Lock Buckets (set/4) ClickSeal biocontainment lids (set/4) 5/7 mL blood tube adapters (set/4, 36 tube capacity) 10 mL blood tube adapters (set/4, 28 tube capacity)



*Indicates a centrifuge, rotor and/or adapter product approved for use in clinical applications; see user manual for details.



Maximize capacity. Higher quality pellets.

There is a rotor available for nearly every application to solve your swinging buckets and fixed angle needs



Swinging buckets

Find your [consumable](#) >

View [rotor details](#) >

Explore [rotor packages](#) >



Fixed angle

Find your [consumable](#) >

View [rotor details](#) >



Swinging buckets rotors



Thermo Scientific Rotor Name	TX-1000	TX-750	BIOShield 1000A	TX-400
Maximum Capacity (mL)	4 x 1000 mL	4 x 750 mL	4 x 250 mL	4 x 400 mL
Maximum Speed (rpm) ¹	4200 rpm	4700 rpm	6000 rpm	5000 rpm
Maximum RCF (x g)	4122 x g	4816 x g	7164 x g	4696 x g
Biocontainment	✓ (Optional)	✓ (Optional)	✓ (Integrated)	✓ (Optional)
Multifuge X4 Pro/Legend XF*	✓	✓	✓	
Megafuge ST4 Plus	✓	✓	✓	
Multifuge X1				✓
Megafuge ST16				✓
Maximum Capacity				
1000 mL Bio-Bottle	4	-	-	-
750 mL Bio-Bottle	4	4	-	-
500 mL Nalgene™ Bottle	4	-	-	-
500 mL Corning™ Bottle	4	4	-	-
250 mL Nalgene Bottle	8	4	4	4
250 mL Corning Bottle	4	4	-	-
225 mL - 175 mL Conical Bottle	8	4	-	4
50 mL Conical Tube	40	28	16	16
50 mL Skirted Conical Tube	20	20	16	-
50 mL Nalgene Oak Ridge Tube	48	28	-	16
30 mL Sterilin™ Universal Tube	28	28	16	12
15 mL Conical Tube	96	56	36	36
10 mL Blood Collection Tube (16 x 100 mm)	148	84	64	56
5/7 mL Blood Collection Tube (13 x 75-100 mm)	196	108	80	76
4.5/6 mL Blood Collection Tube	164	108	-	64
1.5/2 mL Microtube	192	192	224	136
Microplate	24	16	-	-

Biocontainment certification by Public Health England, Porton Down, UK.
¹Maximum speed based on top performance models. Labware not shown to scale.



Thermo Scientific Rotor Name	BIOShield 720	TX-200	HIGHPlate 6000	Fiberlite H3-LV	M-20
Maximum Capacity (mL)	4 x 180 mL	4 x 180 mL	10 Microplates	28 Microplates	6 Microplates
Maximum Speed (rpm) ¹	6300 rpm	5500 rpm	6300 rpm	3600 rpm	4000 rpm
Maximum RCF (x g)	7188 x g	5580 x g	6168 x g	2738 x g	2272 x g
Biocontainment	✓ (Integrated)	✓ (Optional)	✓ (Integrated)		✓ (Optional)
Multifuge X4 Pro/Legend XF*			✓	✓	✓
Megafuge ST4 Plus					✓
Multifuge X1	✓	✓			✓
Megafuge ST16	✓	✓			✓
Maximum Capacity					
1000 mL Bio-Bottle	-	-	-	-	-
750 mL Bio-Bottle	-	-	-	-	-
500 mL Nalgene Bottle	-	-	-	-	-
500 mL Corning Bottle	-	-	-	-	-
250 mL Nalgene Bottle	-	-	-	-	-
250 mL Corning Bottle	-	-	-	-	-
225 mL - 175 mL Conical Bottle	-	-	-	-	-
50 mL Conical Tube	8	4	-	-	-
50 mL Skirted Conical Tube	-	4	-	-	-
50 mL Nalgene Oak Ridge Tube	-	-	-	-	-
30 mL Sterilin Universal Tube	16	4	-	-	-
15 mL Conical Tube	24	20	-	-	-
10 mL Blood Collection Tube (16 x 100 mm)	48	28	-	-	-
5/7 mL Blood Collection Tube (13 x 75-100 mm)	48	32	-	-	-
4.5/6 mL Blood Collection Tube	-	32	-	-	-
1.5/2 mL Microtube	160	48	-	-	-
Microplate	-	-	10	28	6

Biocontainment certification by Public Health England, Porton Down, UK.
¹Maximum speed based on top performance models. Labware not shown to scale.

*Refers to the Legacy model.



Rotors and adapters for swinging buckets rotors



TX-1000 Swinging Buckets Rotor¹, 90°, R_{max} 209 mm

Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4R Pro/Legend XFR*	4200	4122
Multifuge X4 Pro/Legend XF*	3800	3374
MegaFuge ST4R Plus	4200	4122
MegaFuge ST4 Plus	3800	3374
Multifuge X1	n/a	n/a
MegaFuge ST16	n/a	n/a

¹Rotor, buckets and lids sold separately

Adapters



TX-750 Swinging Buckets Rotor¹, 90°, R_{max} 195 mm

Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4 Pro/Legend XFR*	4700	4816
MegaFuge ST4 Plus	4700	4816
Multifuge X1	n/a	n/a
MegaFuge ST16	n/a	n/a

¹Rotor, buckets and lids sold separately

Adapters



Cat. No.	Adapter Image	Description	Rotor Capacity (places x volume, mL)	Maximum Tube Dimensions (Ø x L, mm)
Rotor and Accessories				
75003017		TX-1000 Rotor Cross	4 x 1000 mL or 24 microplate	
75003001		TX-1000 Buckets (Set of 4)		
75007309		TX-1000 ClickSeal Biocontainment Lids (Set of 4)		
75007001		Replacement O-rings		
75007300		1000 mL Bio-Bottle - Polypropylene (Set of 4)		
Adapters for TX-1000 Rotor (Sets of 4)				
75007301		1000 mL Bio-Bottle	4 x 1000 mL	126 x 140
75007304		750 mL Polypropylene Bio-Bottle and any TX-750 Round adapter	4 x 750 mL	98 x 133
75004253		500 mL Nalgene Bottle	4 x 500 mL	70 x 160
75007302		500 mL Corning Bottle (requires Corning #431124)	4 x 500 mL	98 x 133
75005392		250 mL Corning Conical Bottle or 200 mL Nunc™ Bottle or 175 mL Nalgene Conical Bottle	4 x 250 mL	62 x 145
75007305		250 mL Nalgene Bottle	8 x 250 mL	62 x 145
75007305		225 mL BD Falcon™ (requires BD #352090)	8 x 225 mL	62 x 130
75007305		200 mL Nunc Conical Bottle (requires Nunc #377585)	8 x 200 mL	62.5 x 139
75007305		175 mL Nalgene Conical Bottle (requires Nalgene #DS3126-0175)	8 x 175 mL	62 x 130
75003674		50 mL Conical Tube	40 x 50 mL	30 x 120
75004255		50 mL Double Biocontainment Vessel for 50 mL Conical Tube (can be combined with ClickSeal lids)	20 x 50 mL	30 x 120
75004252		50 mL Nalgene Oak Ridge Tube	48 x 50 mL	29.5 x 120
75007306		15 mL Conical Tube	96 x 15 mL	17.5 x 121
75003672		10 mL Blood Collection Tube (16 x 100 mm) or Corex™/Kimble™ Tubes	148 x 10 mL	17 x 110
75003697		9/10 mL Blood Collection Tube (Sarstedt™)	100 x 9/10 mL	17 x 110
75003671		5/7 mL Blood Collection Tube (13 x 75-100 mm)	196 x 5/7 mL	14 x 110
75003709		4.5/6 mL Blood Collection Tube (Greiner™)	164 x 4.5/6 mL	14 x 110
75007303		Microplate Carriers	4 x 6 Standard or 4 x 2 Deep-Well	Height < 88 mm
75003829		Small Bag/Cell Culture Bags 4 x 2 Bags (< 350 mL)	4 x 2 Bags	

Biocontainment certification by Public Health England, Porton Down, UK.

Cat. No.	Adapter Image	Description	Rotor Capacity (places x volume, mL)	Maximum Tube Dimensions (Ø x L, mm)
Rotor and Accessories				
75003180		TX-750 Rotor Cross		
75003608		TX-750 Round Buckets (Set of 4)	4 x 750 mL	
75003609		TX-750 Round ClickSeal Biocontainment Lids (Set of 4)		
75003610		Replacement TX-750 Round O-rings for lids (Set of 4)		
75006443		750 mL Bio-Bottle - Polypropylene (1 each)		
75003795		Microplate Carriers and T-75 Flask (includes tray and pads) (Set of 2)	4 Standard Microplate or 1 Deep well plate per carrier	Maximum height 60 mm
75003617		Microplate Carriers and T-75 Flask (includes tray and pads) (Set of 4)		
Adapters for TX-750 Rotor (Sets of 4, unless stated otherwise)				
Direct		750 mL Bio-Bottle	4 x 750 mL	98 x 133
Direct		500 mL Corning Conical Bottle (requires Corning #431124)	4 x 500 mL	98 x 133
75003792		250 mL Corning Conical Bottle (unsealed buckets only)	4 x 250 mL	62.5 x 139
75003710		250 mL Nalgene Bottle	4 x 250 mL	62.5 x 139
75003710		225 mL/175 mL BD Falcon Conical Bottle (requires BD #352090) (Open buckets only)	4 x 225 mL	62.5 x 139
75003710		200 mL Nunc Conical Bottle (requires Nunc #377585) (Open buckets only)	4 x 200 mL	62.5 x 139
75003710		175 mL Nalgene Conical Bottle (requires Nalgene #DS3126-0175) (Open buckets only)	4 x 175 mL	62.5 x 139
75003713		100 mL Round Bottom Open-Top Tube	8 x 100 mL	45 x 123
75003714		50 mL Conical Tube (includes sealing vessel) can be combined with ClickSeal lids	12 x 50 mL	30 x 120
75003638		50 mL Conical Tube	28 x 50 mL	30 x 120
75003824		50 mL Conical or Skirted Tube	20 x 50 mL	30 x 120
75003715		50 mL Nalgene Oak Ridge Tube	28 x 50 mL	29.5 x 120
75003716		30 mL Sterilin Universal Tube	28 x 30 mL	25.5 x 120
75003639		15 mL Conical Tube	56 x 15 mL	17 x 122
75003719		15 mL Blood Collection Tube (17 x 125 mm) (inner circle only)	32 x 15 mL	17 x 134
75003719		10 mL Blood Collection (16 x 100 mm) or 15 mL Corex/Kimble Tube	84 x 10/15 mL	17 x 117
75003718		14 mL Round or Conical Urine Tube	48 x 14 mL	18 x 124
75003723		5/7 mL or 4.5/6 mL Blood Collection Tube (13 x 75-100 mm)	108 x 5/7 mL	13 x 116
75003724		5 mL RIA or Round Bottom Tube (without cap)	148 x 5 mL	13 x 111
75003732		5/7 mL Round Bottom Tube (without cap) with Decanting Aid	148 x 5/7 mL	13 x 111
75003733		1.5/2 mL Microtube	192 x 2 mL	11 x 50
75008383		T-75 Nunc Easy Flask (Set of 2)	4 x T-75	-
75008384		T-25 Nunc Easy Flask (Set of 2)	8 x T-25	-

Biocontainment certification by Public Health England, Porton Down, UK.

*Refers to the Legacy model.




Rotors and adapters for swinging buckets rotors



TX-750 Swinging Buckets Rotor,
90°, R_{max} 195 mm

Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4R Pro/Legend XFR*, Megafuge ST4R Plus, 230V	4700	4816
Multifuge X4 Pro/Legend XF*, Megafuge ST4 Plus, 230V	4500	4415
Multifuge X4 Pro Series/Legend XF*, Megafuge ST4 Plus Series, 120V	4300	4031
Multifuge X1 or Megafuge ST16	n/a	n/a

Adapters

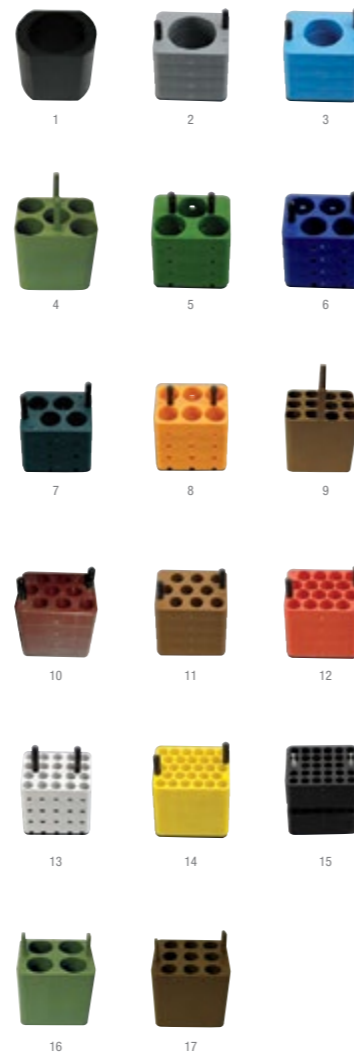

Cat. No.	Adapter Image	Description	Rotor Capacity (places x volume, mL)	Maximum Tube Dimensions (Ø x L, mm)
Rotor and Accessories				
75003180		TX-750 Rotor Cross		
75003614		TX-750 Rectangular Buckets (Set of 4)	4 x 250 mL	
75003615		TX-750 Rectangular Click Seal Biocontainment Lids (Set of 4)		
75003616		Replacement TX-750 Rectangular O-rings for lids (Set of 4)		
Adapters for TX-750 Rotor (Sets of 4)				
75003737	1	250 mL Flat Bottom Bottle	4 x 250 mL	63 x 130
75003738	2	150 mL Round Bottom Open-Top Tube	4 x 150 mL	56 x 124
75003742	3	100 mL Round Bottom Open-Top Tube	4 x 100 mL	46 x 124
75003685	4	50 mL Conical Tube	20 x 50 mL	30 x 121
75003749	5	50 mL Round Bottom Tube	12 x 50 mL	35 x 121
75003750	6	45 mL Flat/Round Tube	16 x 45 mL	30 x 121
75003755	7	30 mL Sterilin Universal Tube	16 x 30 mL	26 x 121
75003756	8	25 mL DIN Round Bottom Tube	24 x 25 mL	26 x 121
75003684	9	15 mL Conical Tube	48 x 15 mL	17 x 120
75003758	10	14 mL Flanged Round Bottom Tube	36 x 14 mL	18 x 121
75003759	11	14 mL Round or Conical Urine Tube	36 x 15 mL	17 x 124
75003767	12	10 mL Blood Collection (16 x 100 mm) or 15 mL DIN Tube	64 x 10/15 mL	17 x 118
75003768	13	5/7 mL or 4.5/6 mL Blood Collection Tube	80 x 5/7 mL	14 x 118
75003769	14	5/7 mL Round Bottom Tube (13 x 75-100 mm)	112 x 5/7 mL	13 x 119
75003770	15	1.5/2 mL Microtube	224 x 2 mL	11 x 50



BIOShield 1000A Swinging Buckets Rotor, 90°, R_{max} 178 mm

Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4 Pro/Legend XF*	6000	7164
Megafuge ST4 Plus	5300	5590
Multifuge X1	n/a	n/a
Megafuge ST16	n/a	n/a

Cat. No.	Adapter Image	Description	Rotor Capacity (places x volume, mL)	Maximum Tube Dimensions (Ø x L, mm)
Rotor and Accessories				
75003182		BIOShield 1000A	4 x 250 mL	
Adapters for BIOShield 1000A Rotor (Sets of 4)				
75003737	1	250 mL Flat Bottom Bottle	4 x 250 mL	63 x 130
75003738	2	150 mL Round Bottom Open-Top Tube	4 x 150 mL	56 x 124
75003742	3	100 mL Round Bottom Open-Top Tube	4 x 100 mL	46 x 124
75003749	5	50 mL Round Bottom Tube	12 x 50 mL	35 x 121
75003750	6	45 mL Flat/Round Tube	16 x 45 mL	30 x 121
75003755	7	30 mL Sterilin Universal Tube	16 x 30 mL	30 x 121
75003756	8	25 mL DIN Round Bottom Tube	24 x 25 mL	26 x 121
75003758	10	14 mL Flanged Round Bottom Tube	36 x 14 mL	18 x 121
75003759	11	14 mL Round or Conical Urine Tube	36 x 15 mL	17 x 124
75003767	12	10 mL Blood Collection (16 x 100 mm) or 15 mL DIN Tube	64 x 10/15 mL	17 x 118
75003768	13	5/7 mL or 4.5/6 mL Blood Collection Tube	80 x 5/7 mL	14 x 118
75003769	14	5/7 mL Round Bottom Tube (13 x 75-100 mm)	112 x 5/7 mL	13 x 119
75003770	15	1.5/2 mL Microtube	224 x 2 mL	11 x 50
75003643	16	50 mL Conical Tube	12 x 50 mL	30 x 121
75003642	17	15 mL Conical Tube	36 x 15 mL	17 x 120

Biocontainment certification by Public Health England, Porton Down, UK.



TX-400 Swinging Buckets Rotor¹,
90°, R_{max} 168 mm

Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4 Pro/Legend XF*	n/a	n/a
Megafuge ST4 Plus	n/a	n/a
Multifuge X1	5000	4696
Megafuge ST16	5000	4696

Adapters



Cat. No.	Adapter Image	Description	Rotor Capacity (places x volume, mL)	Maximum Tube Dimensions (Ø x L, mm)
Rotor and Accessories				
75003181		TX-400 Rotor Cross		
75003655		TX-400 Round Buckets (Set of 4)	4 x 400 mL	
75003656		TX-400 Round ClickSeal Biocontainment Lids (Set of 4)		
75003657		Replacement TX-400 O-rings for lids (Set of 4)		
75007585		400 mL Bio-Bottle - Polypropylene (Set of 12)		
Adapters for TX-400 Rotor (Sets of 4)				
Direct		400 mL Bio-Bottle	4 x 400 mL	80 x 125
75003788	1	250 mL Thermo Scientific Nalgene™ Bottle	4 x 250 mL	62 x 135
75003788	1	200 mL Thermo Scientific Nunc Conical Bottle (requires Nunc #377585)	4 x 200 mL	62 x 125
75003788	1	225 mL/175 mL BD Falcon Conical Bottle (requires BD #352090)	4 x 175 mL	62 x 125
75003708	2	100 mL Round Bottom Open-Top Tube	4 x 100 mL	45 x 117
75003707	3	50 mL DIN Round Bottom Tube	12 x 50 mL	34.5 x 105
75003799	4	50 mL Nalgene™ Oak Ridge Tube	16 x 50 mL	29.5 x 114
75003683	5	50 mL Conical Tube	16 x 50 mL	30 x 116
75003703	6	30/25 mL DIN Round/Flat Bottom Tube	20 x 30 mL	25.5 x 108
75003706	7	30 mL Sterilin Universal Tube	12 x 30 mL	25 x 110
75003682	8	15 mL Conical Tube	36 x 15 mL	17 x 121
75003794	9	15 mL Blood Collection Tube (17 x 125 mm)	16 x 15 mL	17 x 131
75003704	10	15 mL Round Bottom Tube (Sarstedt)	40 x 15 mL	17 x 105
75003798	11	14 mL Round or Conical Urine Tube	28 x 14 mL	17 x 112
75003681	12	10 mL Blood Collection (16 x 100 mm)	56 x 10/15 mL	16 x 113
75003680	13	5/7 mL Blood Collection Tube (13 x 75-100 mm)	76 x 5/7 mL	13 x 110
75003825	14	4.5/6 mL Blood Collection Tube (Greiner)	64 x 4.5/6 mL	13 x 110
75003793	15	3 mL RIA or Round Bottom Tube (without cap)	76 x 3/5 mL	11 x 110
75003700	16	1.5/2 mL Microtube	136 x 2 mL	11 x 45

Biocontainment certification by Public Health England, Porton Down, UK.

more >

¹Refers to the Legacy model.



Rotors and adapters for swinging buckets rotors



BIOShield 720 Swinging Buckets Rotor, 90°, R_{max} 162 mm

Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4 Pro/Legend XF*	n/a	n/a
Megaforce ST4 Plus	n/a	n/a
Multifuge X1	6300	7188
Megaforce ST16	5300	5088

Cat. No.	Adapter Image	Description	Rotor Capacity (places x volume, mL)	Maximum Tube Dimensions (Ø x L, mm)
Rotor and Accessories				
75003183		BIOShield 720	4 x 180 mL	n/a
75003693		180 mL Bio-Bottle - Polypropylene (Set of 12)		
BIOShield 720A Adapters (Sets of 4)				
75003813		150 mL Round Bottom Open-Top Tube	4 x 150 mL	50 x 122
75003814		100 mL Round Bottom Open-Top Tube	4 x 100 mL	44 x 120
75003677		50 mL Conical Tube	8 x 50 mL	30 x 116
75003816		50 mL DIN Round Bottom Tube	8 x 50 mL	34 x 109
75003818		30 mL Sterilin Universal Tube	16 x 30 mL	26 x 105
75003817		25 mL DIN Round Bottom Tube	20 x 25 mL	25 x 107
75003678		15 mL Conical Tube	24 x 15 mL	17 x 120
75003820		15 mL Sarstedt Tube	44 x 15 mL	17 x 103
75003701		10 mL Blood Collection Tube (16 x 100 mm)	48 x 10 mL	17 x 105
75003821		5/7 mL Blood Collection Tube (13 x 75-100 mm)	48 x 5/7 mL	13 x 106
75003822		5/7 mL Round Bottom Open-Top Tube	80 x 5/7 mL	13 x 104
75003823		1.5/2 mL Conical Microtube	160 x 2 mL	11 x 45

Biocontainment certification by Public Health England, Porton Down, UK.



TX-200 Swinging Buckets Rotor¹, 90°, R_{max} 165 mm

Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4 Pro/Legend XF*	n/a	n/a
Megaforce ST4 Plus	n/a	n/a
Multifuge X1 R 230V and 120V	5500	5580
Multifuge X1 120V	5250	5084
Megaforce ST16 R 230V and 120V	5500	5580
Megaforce ST16 120V	5250	5084

¹Rotor, buckets and lids sold separately

Cat. No.	Adapter Image	Description	Rotor Capacity (places x volume, mL)	Maximum Tube Dimensions (Ø x L, mm)
Rotor and Accessories				
75003658		TX-200 Rotor Cross		
75003659		TX-200 Round Buckets (Set of 4)	4 x 180 mL	
75003660		TX-200 Round ClickSeal Biocontainment Lids (Set of 4)		
75003687		Replacement TX-200 O-rings for lids (Set of 4)		
75003800		180 mL Bio-Bottle - Polypropylene (Set of 12)		
Adapters for TX-200 Rotor (Sets of 4)				
Direct		180 mL Bio-Bottle	4 x 180 mL	56.5 x 123
75003801		100 mL Round Bottom Open-Top Tube	4 x 100 mL	48 x 122
75003802		50 mL DIN Round Bottom Tube	4 x 50 mL	35 x 126
75003815		50 mL Round Bottom Tube	4 x 50 mL	32 x 126
75003803		50 mL Conical or Skirted Tube	4 x 50 mL	30 x 126
75003804		30 mL Sterilin™ Universal Tube	4 x 30 mL	26 x 118
75003805		25 mL DIN Round/Flat Bottom Tube	8 x 25 mL	25 x 118
75003806		20 mL Round Bottom Tube	12 x 25 mL	22 x 125
75003771		15 mL Conical Tube	20 x 20 mL	17 x 121
75003809		15 mL Blood Collection Tube	16 x 15 mL	17 x 121
75003807		14 mL Conical Urine Tube	8 x 14 mL	17 x 114
75003808		10 mL Blood Collection (16 x 100 mm) or 15 mL Cores/Kimble Tube	28 x 10/15 mL	17 x 121
75003785		5/7 mL or 4.5/6 mL Blood Collection Tube (13 x 75-100 mm)	32 x 5/7 mL	13 x 122
75003810		5/7 mL Round Bottom Open-Top Tube	36 x 5/7 mL	13 x 118
75003811		3/5 mL RIA or Round Bottom Tube (without cap)	48 x 3/5 mL	12 x 124
75003812		1.5/2 mL Microtube	48 x 2 mL	11 x 45

Biocontainment certification by Public Health England, Porton Down, UK.



more >

*Refers to the Legacy model.



Rotors and adapters for swinging buckets rotors

page 4 of 4

**HIGHPlate 6000 Swinging Buckets Rotor, 90°, R_{max} 139 mm**

Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4 Pro/Legend XF*	6300	6168
Megafuge ST4 Plus	n/a	n/a
Multifuge X1	n/a	n/a
Megafuge ST16	n/a	n/a

Cat. No.	Description	Rotor Capacity (places x volume, mL)	Maximum Height
Rotor			
75003606	HIGHPlate 6000 Rotor	2 x 5 Standard Microplate or 2 x 1 Filter plates	< 80 mm

**Fiberlite H3-LV Swinging Buckets Rotor, 90°, R_{max} 189 mm**

Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4 Pro/Legend XF*	3600	2738
Megafuge ST4 Plus	n/a	n/a
Multifuge X1	n/a	n/a
Megafuge ST16	n/a	n/a

Cat. No.	Description	Rotor Capacity (places x volume, mL)	Maximum Height
Rotor			
75003665	Fiberlite H3-LV Rotor	4 x 7 Standard Microplate or 2 x 4 Deepwell	< 100 mm

**M-20 Swinging Buckets Rotor¹, 90°, R_{max} 127 mm**

Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4 Pro/Legend XF*	4000	2272
Megafuge ST4 Plus	4000	2272
Multifuge X1	4000	2272
Megafuge ST16	4000	2272

¹ Biocontainment system sold separately

Cat. No.	Description	Rotor Capacity (places x volume, mL)	Maximum Height
Rotor and Accessories			
75003624	M-20 Microplate Rotor including standard open carriers	2 x 3 Standard Microplate	< 65 mm
75003625	ClickSeal Carriers including Carrier and lids	2 x 1 Deepwell	< 45 mm
75002011	ClickSeal Replacement lids (Set of 2)		
75002012	Replacement M-20 O-rings for lids		

Biocontainment certification by Public Health England, Porton Down, UK.

*Refers to the Legacy model.



Fixed angle rotors



Thermo Scientific Rotor Name	Fiberlite F14-6x250 LE	Fiberlite F15-6x100y	Fiberlite F10-6x100 LEX	Fiberlite F13-14x50cy	Fiberlite F15-8x50cy	CLINIConic
Maximum Capacity (mL)	6 x 250 mL	6 x 100 mL	6 x 100 mL	14 x 50 mL	8 x 50 mL	30 x 15 mL
Maximum Speed (rpm) ¹	11000 rpm	15000 rpm	10000 rpm	10500 rpm	14500 rpm	5650 rpm
Maximum RCF (x g)	18533 x g	24652 x g	17105 x g	15038 x g	24446 x g	4997 x g
Biocontainment	✓	✓	✓	✓	✓	✓
Multifuge X4 Pro/Legend XF*	✓	✓	✓	✓	✓	✓
Megafuge ST4 Plus	✓	✓	✓	✓	✓	✓
Multifuge X1	✓	✓	✓	✓	✓	✓
Megafuge ST16	✓	✓	✓	✓	✓	✓
Maximum Capacity	250 mL Nalgene Bottle	6	-	-	-	-
	85 mL Nalgene Oak Ridge Tube	6	6	6	-	-
	50 mL Nalgene Oak Ridge Tube	6	6	6	14	8
	50 mL Conical Tube	6	6	6	14	8
	30 mL Nalgene Oak Ridge Tube	12	6	6	14	8
	16 mL Nalgene Oak Ridge Tube	30	6	6	14	8
	15 mL Conical Tube	30	6	6	14	8
	10 mL Nalgene Oak Ridge Tube	42	12	12	14	8
	10 mL Blood Collection Tube (16 x 100 mm)	-	-	-	-	30
	5/7 mL Blood Collection Tube (13 x 75-100 mm)	-	-	-	-	30
	6.5 mL Round Tube	-	12	12	-	-
	1.5/2 mL Microtube	-	24	24	-	-

Biocontainment certification by Public Health England, Porton Down, UK.
¹Maximum speed based on top performance models. Labware not shown to scale.



Thermo Scientific Rotor Name	8 x 50 Sealed	HIGHConic II	MicroClick 18 x 5	Fiberlite F21-48x2	Microliter 48 x 2	MicroClick 30 x 2	Microliter 30 x 2
Maximum Capacity (mL)	8 x 50 mL	6 x 100 mL	18 x 5 mL	48 x 2 mL	48 x 2 mL	30 x 2 mL	30 x 2 mL
Maximum Speed (rpm) ¹	6700 rpm	10350 rpm	15000 rpm	15200 rpm	15200 rpm	14000 rpm	15200 rpm
Maximum RCF (x g)	7177 x g	15090 x g	24652 x g	25055 x g	25314 x g	21693 x g	25830 x g
Biocontainment	✓ (Individual)	✓	✓	✓	✓	✓	✓
Multifuge X4 Pro/Legend XF*	✓	✓	✓	✓	✓	✓	✓
Megafuge ST4 Plus	✓	✓	✓	✓	✓	✓	✓
Multifuge X1	✓	✓	✓	✓	✓	✓	✓
Megafuge ST16	✓	✓	✓	✓	✓	✓	✓
Maximum Capacity	250 mL Nalgene Bottle	-	-	-	-	-	-
	85 mL Nalgene Oak Ridge Tube	-	6	-	-	-	-
	50 mL Nalgene Oak Ridge Tube	-	6	-	-	-	-
	50 mL Conical Tube	8	6	-	-	-	-
	30 mL Nalgene Oak Ridge Tube	-	6	-	-	-	-
	16 mL Nalgene Oak Ridge Tube	-	6	-	-	-	-
	15 mL Conical Tube	8	6	-	-	-	-
	10 mL Nalgene Oak Ridge Tube	-	12	-	-	-	-
	10 mL Blood Collection Tube (16 x 100 mm)	8	-	-	-	-	-
	5/7 mL Blood Collection Tube (13 x 75-100 mm)	8	-	-	-	-	-
	6.5 mL Round Tube	-	12	-	-	-	-
	5 mL microtube/conical	-	-	18	-	-	-
	1.5/2 mL Microtube	-	24	18	48	48	30 (Including spin columns)

Biocontainment certification by Public Health England, Porton Down, UK.
¹Maximum speed based on top performance models. Labware not shown to scale.

*Refers to the Legacy model.



Rotors and adapters for fixed angle rotors



Fiberlite F14-6x250 LE Fixed Angle Rotor, 23°, R_{max} 137 mm

Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4 Pro/Legend XF* 230V	11000	18533
Multifuge X4 Pro/Legend XF* 120V	10000	15317
Megafuge ST4 Plus	n/a	n/a
Multifuge X1	n/a	n/a
Megafuge ST16	n/a	n/a

Cat. No.	Description	Rotor Capacity (places x volume, mL)	Maximum Tube Dimensions (Ø x L, mm)
Rotor and Accessories			
75003662	Fiberlite F14-6 x 250LE Rotor for 250 mL Nalgene Bottles	6 x 250 mL	62 x 125
021-062153	Replacement O-ring Kit		
3141-0250	250 mL Nalgene Polypropylene Bottle (Set of 36)		
3140-0251	250 mL Nalgene Polycarbonate Bottle (Set of 36)		
Adapters for F14-6 X 250LE Rotor (Sets of 2, unless stated otherwise)			
010-1119	85 mL Nalgene Oak Ridge Tube	6 x 85 mL	38 x 109
010-0138	50 mL Nalgene Oak Ridge Tube	6 x 50 mL	29.5 x 121
75100136	50 mL Conical Tube	6 x 50 mL	30 x 121
010-1072	30 mL Nalgene Oak Ridge Tube / 38 mL Round Tube	12 x 30 mL	26 x 105
010-1074	16 mL Nalgene Oak Ridge Tube	30 x 16 mL	18 x 107
010-1410	15 mL Conical Tube	30 x 15 mL	17 x 121
00389	10 mL Nalgene Oak Ridge Tube / 12 mL Round Tube (each)	42 x 10 mL	17 x 110



Fiberlite F10-6x100 LEX Fixed Angle Rotor, 45°, R_{max} 122 mm

Fiberlite F15-6x100y Fixed Angle Rotor, 25°, R_{max} 98 mm

Centrifuge	F10 Maximum Speed (rpm)	F10 Maximum RCF (x g)	F15 Maximum Speed (rpm)	F15 Maximum RCF (x g)
Multifuge X4 Pro/Legend XF*	10500	15038	15000	24652
Megafuge ST4 Plus	10500	15038	13000	18516
Multifuge X1	10500	15038	15000	24652
Megafuge ST16	10500	15038	13000	18516

Cat. No.	Description	Rotor Capacity (places x volume, mL)	Maximum Tube Dimensions (Ø x L, mm)
Rotor and Accessories			
75003698	Fiberlite F15-6 x 100y Rotor for 100 mL/85 mL Nalgene	6 x 100 mL	38 x 106
75003340	Fiberlite F10-6 x 100 LEX Rotor for 100mL/85mL Nalgene	6 x 100 mL	38 x 106
021-069031	Replacement O-ring Kit		
3118-0085	85 mL Nalgene Oak Ridge Tube (Set of 50)	6 x 85 mL	38 x 106
Adapters for F15-6 x 100y Rotor (Each, unless stated otherwise)			
75003102	50 mL Nalgene Oak Ridge Tube	6 x 50 mL	29.5 x 117
75003103	50 mL Conical Tube	6 x 50 mL	30 x 121
75003094	30 mL Nalgene™ Oak Ridge Tube / 38 mL Round Tube	6 x 30 mL	26 x 110
76002906	16 mL Nalgene Oak Ridge Tube	6 x 16 mL	18 x 123
75003095	15 mL Conical Tube	6 x 15 mL	17 x 123
75003093	10 mL Nalgene Oak Ridge Tube / 12 mL Round Tube	12 x 12 mL	17 x 95
75003092	6.5 mL Round Bottom Tube	12 x 6.5 mL	13 x 114
75003091	1.5/2 mL Microtube	24 x 2 mL	11 x 40

Biocontainment certification by Public Health England, Porton Down, UK.



Fiberlite F13-14x50cy Fixed Angle Rotor, 34°, R_{max} 153 mm

Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4 Pro/Legend XF* 230V	10000	17105
Multifuge X4 Pro/Legend XF* 120V	9250	14636
Megafuge ST4 Plus	n/a	n/a
Multifuge X1	8500	12359
Megafuge ST16	n/a	n/a

Cat. No.	Description	Rotor Capacity (places x volume, mL)	Maximum Tube Dimensions (Ø x L, mm)
Rotor and Accessories			
75003661	Fiberlite F13-14x50cy for 50 mL Conical Tubes	14 x 50 mL	30 x 121
021-149027	Replacement O-ring Kit		
Adapters for F13-14x50cy Rotor (Set of 2, unless otherwise stated)			
010-0377	50 mL Nalgene Oak Ridge Tube	14 x 50 mL	29.5 x 108
010-1147	30 mL Nalgene Oak Ridge Tube	14 x 30 mL	26 x 102
010-0376	16 mL Nalgene Oak Ridge Tube	14 x 16 mL	18 x 123
75100378	15 mL Conical Tube	14 x 15 mL	17 x 121
010-1340	15 mL Nunc Conical Tube and Filtration Devices	14 x 15 mL	17 x 121
010-1311	10 mL Nalgene Oak Ridge Tube / 12 mL Round Tube	28 x 12 mL	17 x 95



Fiberlite F15-8x50cy Fixed Angle Rotor, 25°, R_{max} 104 mm

Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4 Pro/Legend XF*	14500	24446
Megafuge ST4 Plus	n/a	n/a
Multifuge X1 230V	14500	24446
Multifuge X1 120V	14000	22789
Megafuge ST16	n/a	n/a

Cat. No.	Description	Rotor Capacity (places x volume, mL)	Maximum Tube Dimensions (Ø x L, mm)
Rotor and Accessories			
75003663	Fiberlite F15-8 x 50cy for 50 mL Conical Tubes	8 x 50 mL	30 x 121
021-085077	Replacement O-ring Kit		
Adapters for F15-8x50cy Rotor (Set of 2, unless otherwise stated)			
010-0377	50 mL Nalgene Oak Ridge Tube	8 x 50 mL	29 x 108
010-1147	30 mL Nalgene Oak Ridge Tube	8 x 30 mL	26 x 102
010-0376	16 mL Nalgene Oak Ridge Tube	8 x 16 mL	18 x 123
75100378	15 mL Conical Tube	8 x 15 mL	17 x 121
010-1340	15 mL Nunc Conical Tube and Filtration Devices	8 x 15 mL	17 x 121
010-1311	10 mL Nalgene Oak Ridge Tube / 12 mL Round Tube	8 x 12 mL	17 x 95

Biocontainment certification by Public Health England, Porton Down, UK.

*Refers to the Legacy model.



Rotors and adapters for fixed angle rotors



CLINIConic Fixed Angle Rotor, 37°, R_{max} 140 mm

Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4 Pro/Legend XF*	n/a	n/a
Megafuge ST4 Plus	n/a	n/a
Multifuge X1	5650	4997
Megafuge ST16	5650	4997

Cat. No.	Description	Rotor Capacity (places x volume, mL)	Maximum Tube Dimensions (Ø x L, mm)
Rotor			
75003623	CLINIConic Fixed Angle Rotor for 15 mL Round and Conical Tubes	30 x 15 mL	17 x 121
Adapters for CLINIConic Rotor (Each)			
75003702	10 mL Round Bottom Tube	30 x 10 mL	17 x 131
11172596	7 mL Blood Collection Tubes (13 x 100 mm)	30 x 7 mL	17 x 131
11172595	5 mL Blood Collection Tubes (13 x 75 mm)	30 x 5 mL	13 x 95



8 x 50 Individually Sealed Fixed Angle Rotor, 45°, R_{max} 143 mm

Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4 Pro/Legend XF*	n/a	n/a
Megafuge ST4 Plus	n/a	n/a
Multifuge X1 230V	6700	7177
Multifuge X1 120V	6500	6755
Megafuge ST16 230V	6700	7177
Megafuge ST16 120V	6500	6755

Cat. No.	Description	Rotor Capacity (places x volume, mL)	Maximum Tube Dimensions (Ø x L, mm)
Rotor and Accessories			
75003694	8 x 50 Individually Sealed for 50 mL Conical Tubes	8 x 50 mL	30 x 121
75003689	Replacement O-ring Kit		
Adapters for 8 x 50 Sealed Rotor (Set of 6)			
75005755	15 mL Conical Tube	8 x 50 mL	17 x 121
75005747	10 mL Blood Collection Tube (16 x 100 mm)	8 x 10 mL	16 x 110
75005748	7 mL Blood Collection Tube (13 x 100 mm)	8 x 7 mL	13 x 110
75005749	3.5 mL Blood Collection Tube	16 x 3.5 mL	11 x 75

Biocontainment certification by Public Health England, Porton Down, UK.



HIGHConic II Aluminium Fixed Angle Rotor, 45°, R_{max} 126 mm

Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4 Pro/Legend XF*	10350	15090
Megafuge ST4 Plus	10350	15090
Multifuge X1	10350	15090
Megafuge ST16	10350	15090

Cat. No.	Description	Rotor Capacity (places x volume, mL)	Maximum Tube Dimensions (Ø x L, mm)
Rotor and Accessories			
75003620	HIGHConic Fixed Angle Rotor for 100 mL/85 mL Nalgene Oak Ridge Tube	6 x 100 mL	38 x 123
75003058	Replacement O-ring Kit		
3118-0085	85 mL Nalgene Oak Ridge Tube (Set of 50)	6 x 85 mL	38 x 106
Adapters for HIGHConic Rotor (Each, unless stated)			
75003102	50 mL Nalgene Oak Ridge Tube	6 x 50 mL	30 x 117
75003103	50 mL Conical Tube	6 x 50 mL	30 x 121
75003094	30 mL Nalgene Oak Ridge Tube / 38 mL Round Tube	6 x 30 mL	26 x 110
76002906	16 mL Nalgene Oak Ridge Tube	6 x 16 mL	18 x 123
75003095	15 mL Conical Tube	6 x 15 mL	17 x 123
75003093	10 mL Nalgene Oak Ridge Tube / 12 mL Round Tube	12 x 12 mL	17 x 95
75003092	6.5 mL Round Bottom Tube	12 x 6.5 mL	13 x 114
75003091	1.5/2 mL Microtube	24 x 2 mL	11 x 40



MicroClick 18 x 5 Aluminum Fixed Angle Rotor, 45°, R_{max} 102 mm

Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4 Pro/Legend XF*	15000	24652
Megafuge ST4 Plus	15000	24652
Multifuge X1	15000	24652
Megafuge ST16	15000	24652

Cat. No.	Description	Rotor Capacity (places x volume, mL)	Maximum Tube Dimensions (Ø x L, mm)
Rotor and Accessories			
75005765	18 x 5 MicroClick Rotor for 5mL Conical and Flip Top Tubes	18 x 5 mL	17 x 6
20059119	ClickSeal	n/a	n/a
75005726	O-Ring Set	n/a	n/a
75003500	O-Ring Grease	n/a	n/a
Adapters for MicroClick Rotor			
75005756	1.5/2 mL Microtube	18 x 5	11 x 45



Fiberlite F21-48 x 2 Fixed Angle Rotor, 45°, R_{max} 97 mm

Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4 Pro/Legend XF*	15200	25055
Megafuge ST4 Plus	15200	25055
Multifuge X1	15200	25055
Megafuge ST16	15200	25055

Cat. No.	Description	Rotor Capacity (places x volume, mL)	Maximum Tube Dimensions (Ø x L, mm)
Rotor and Accessories			
75003664	Fiberlite F21-48 x 2 for 1.5/2 mL Microtubes	48 x 2 mL	11 x 50
021-489021	Replacement O-ring Kit		

Biocontainment certification by Public Health England, Porton Down, UK.

*Refers to the Legacy model.



Rotors and adapters for fixed angle rotors

page 3 of 3

**Microliter 48 x 2 Aluminum Fixed Angle Rotor, 45°, R_{max} 98 mm**

Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4 Pro/Legend XF*	15200	25314
Megafuge ST4 Plus	15200	25314
Multifuge X1	15200	25314
Megafuge ST16	15200	25314

Cat. No.	Description	Rotor Capacity (places x volume, mL)	Maximum Tube Dimensions (Ø x L, mm)
Rotor and Accessories			
75003602	48 x 2 Aluminum for 1.5/2 mL Microtubes	48 x 2 mL	11 x 50
75003349	Replacement O-ring Kit		

**MicroClick 30 x 2 Aluminum Fixed Angle Rotor, 45°, R_{max} 99 mm**

Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4 Pro/Legend XF*	14000	21693
Megafuge ST4 Plus	14000	21693
Multifuge X1	14000	21693
Megafuge ST16	14000	21693

Cat. No.	Description	Rotor Capacity (places x volume, mL)	Maximum Tube Dimensions (Ø x L, mm)
Rotor and Accessories			
75005719	30 x 2 MicroClick for 1.5/2 mL Microtubes and Spin Columns	30 x 2 mL	11 x 50
75005726	Replacement O-ring Kit		

**Microliter 30 x 2 Aluminum Fixed Angle Rotor, 45°, R_{max} 100 mm**

Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4 Pro/Legend XF*	15200	25830
Megafuge ST4 Plus	15200	25830
Multifuge X1	15200	25830
Megafuge ST16	15200	25830

Cat. No.	Description	Rotor Capacity (places x volume, mL)	Maximum Tube Dimensions (Ø x L, mm)
Rotor and Accessories			
75003652	30 x 2 Aluminum for 1.5/2 mL Microtubes	30 x 2 mL	11 x 50
75003349	Replacement O-ring Kit		

Biocontainment certification by Public Health England, Porton Down, UK.

Microliter Rotor Adapters

Cat. No.	Description	Maximum Tube Dimensions (Ø x L, mm)
Adapters for all Microliter Rotors (Set of 30)		
75005753	0.5 mL Microtube	8 x 44
75005754	0.25 mL Microtube	6 x 46
75005752	0.2 mL PCR Tube	5 x 20

*Refers to the Legacy model.



Flexible combinations

Whether your work is in healthcare, research or industrial applications, find the flexible combination you need for both basic and advanced centrifuge performance levels

[Go to research applications](#)[Go to clinical applications*](#)

*Please refer to centrifuge, rotor and adapter product details for products approved for use in clinical applications.



Cell and tissue culture applications

- Conduct multiple research protocols on one centrifuge
- Obtain faster results in our highest capacity TX-1000 rotor
- Compliant with the latest global regulatory and safety standards

[View full cell culture packages](#)

[View rotor-only cell culture packages](#)





Clinical applications*

- **Faster separations** with swinging buckets rotor speeds up to 6300 rpm / 7188 x g
- Process more samples with the ability to load **196 blood tubes per run**
- **Auto-Lock™** rotor exchange simplifies cleaning to **limit potential biohazard** contamination
- Programming options include up to 100 saved programs for **repeatability** by even the newest laboratory members and **two levels of security**, ideal for controlled and regulated environments
- **Enables compliance** to the latest global clinical and **safety standards**

[View full clinical packages](#)

[View rotor-only clinical packages](#)

*Please refer to centrifuge, rotor and adapter product details for products approved for use in clinical applications.



© 2019 Thermo Fisher Scientific Inc. All rights reserved. All trademarks are the property of Thermo Fisher Scientific and its subsidiaries unless otherwise specified.

BRGPPROMULTMEGA 0719



Ready to print?

Click the button below to print this entire brochure, along with its supplement catalog pages.

Only looking to print the supplemental catalog? Simply choose the print page range option, and select pages 41–56.

[Print catalog](#)

Rotors and adapters for 90° Rotor




Cat. No.	Adapter Image	Description
75003180		TX-750 Rotor Cross
75003808		TX-750 Round Buckets (Set of 4)
75003809		TX-750 Round ClickSeal Biocontainment Lids (Set of 4)
75003810		Replacement TX-750 Round O-rings for lids (Set of 4)
75006443		750 mL Bio-Bottle - Polypropylene (1 each)
75003795		Microplate Carriers and T-75 Flask (includes 2 pads) (Set of 2)
75003817		Microplate Carriers and T-75 Flask (includes 2 pads) (Set of 4)
Adapters for TX-750 Rotor (Set of 4, unless stated otherwise)		
Direct		750 mL Bio-Bottle
75003792		250 mL Corning Conical Bottle (requires Corning bucket)
75003710		250 mL Nalgene Bottle
75003710		225 mL/175 mL BD Falcon Conical Bottle (requires BD #552090) (Open buckets only)
75003710		250 mL Nunc Conical Bottle (requires Nunc bucket) (Open buckets only)
75003710		175 mL Nalgene Conical Bottle (requires Nalgene #DS126-0175) (Open buckets only)
75003713		100 mL Round Bottom Open-Top Tube
75003714		50 mL Conical Tube (includes sealing vessel can be combined with ClickSeal lids)
75003838		50 mL Conical Tube
75003824		50 mL Conical or Skirted Tube
75003715		50 mL Nalgene Oak Ridge Tube
75003716		30 mL Sterilin Universal Tube
75003839		15 mL Conical Tube
75003719		15 mL Blood Collection Tube (17 x 125 mm inner circle only)
75003719		15 mL Blood Collection (16 x 100 mm) or 15 mL Corning/Kimble Tube
75003718		14 mL Round or Conical Urine Tube
75003723		5/7 mL, or 4.5/6 mL Blood Collection Tube
75003724		5 mL RIA or Round Bottom Tube (without Decanting Aid)
75003732		5/7 mL Round Bottom Tube (without Decanting Aid)
75003733		1.5/2 mL Microtube
75003833		T-75 Nunc Easy Flask (Set of 2)
75003834		T-25 Nunc Easy Flask (Set of 2)

*Refers to the Legacy model.

Thermo Scientific™ MULTIFUGE™ X PRO AND MEGAFUGE™ ST PLUS CENTRIFUGE SERIES

Swinging buckets rotors



Thermo Scientific Rotor Name	BIOShield 720	TX-200	HIGHPlate 6000	Fiberlite H3-LV	M-20
Maximum Capacity (mL)	4 x 180 mL	4 x 180 mL	10 Microplates	28 Microplates	6 Microplates
Maximum Speed (rpm)	6300 rpm	5500 rpm	6300 rpm	3600 rpm	4000 rpm
Maximum RCF (x g)	7188 x g	5580 x g	6168 x g	2738 x g	2272 x g
Biocontainment	✓ (integrated)	✓ (optional)	✓ (integrated)	✓ (optional)	✓ (optional)
Multifuge X4 Pro Legend XP*	✓	✓	✓	✓	✓
Multifuge ST4 Plus	✓	✓	✓	✓	✓
Multifuge X1	✓	✓	✓	✓	✓
Multifuge ST16	✓	✓	✓	✓	✓
Maximum Capacity	1000 mL Bio-Bottle	750 mL Bio-Bottle	500 mL Nalgene Bottle	500 mL Corning Bottle	250 mL Nalgene Bottle
	250 mL Corning Bottle	225 mL - 175 mL Conical Bottle	50 mL Conical Tube	50 mL Skirted Conical Tube	50 mL Nalgene Oak Ridge Tube
	30 mL Sterilin Universal Tube	15 mL Conical Tube	10 mL Blood Collection Tube (16 x 100 mm)	5/7 mL Blood Collection Tube (13 x 75-100 mm)	4.5/6 mL Blood Collection Tube
	1.5/2 mL Microtube	160	48	32	48
Microplate	10	28	6		

*Refers to the Legacy model.

Thermo Scientific™ MULTIFUGE™ X PRO AND MEGAFUGE™ ST PLUS CENTRIFUGE SERIES

ThermoFisher
SCIENTIFIC

Swinging buckets rotors



Thermo Scientific Rotor Name	TX-1000	TX-750	BIOShield 1000A	TX-400
Maximum Capacity (mL)	4 x 1000 mL	4 x 750 mL	4 x 250 mL	4 x 400 mL
Maximum Speed (rpm) ¹	4200 rpm	4700 rpm	6000 rpm	5000 rpm
Maximum RCF (x g)	4122 x g	4816 x g	7164 x g	4696 x g
Biocontainment	✓ (Optional)	✓ (Optional)	✓ (Integrated)	✓ (Optional)
Multifuge X4 Pro/Legend XF*	✓	✓	✓	✓
Megafuge ST4 Plus	✓	✓	✓	✓
Multifuge X1				✓
Megafuge ST16				✓
Maximum Capacity	1000 mL Bio-Bottle	4	–	–
	750 mL Bio-Bottle	4	4	–
	500 mL Nalgene™ Bottle	4	–	–
	500 mL Corning™ Bottle	4	4	–
	250 mL Nalgene Bottle	8	4	4
	250 mL Corning Bottle	4	4	–
	225 mL - 175 mL Conical Bottle	8	4	4
	50 mL Conical Tube	40	28	16
	50 mL Skirted Conical Tube	20	20	16
	50 mL Nalgene Oak Ridge Tube	48	28	16
	30 mL Sterilin™ Universal Tube	28	28	16
	15 mL Conical Tube	96	56	36
	10 mL Blood Collection Tube (16 x 100 mm)	148	84	64
	5/7 mL Blood Collection Tube (13 x 75-100 mm)	196	108	80
	4.5/6 mL Blood Collection Tube	164	108	64
	1.5/2 mL Microtube	192	192	224
	Microplate	24	16	–














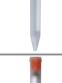
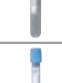





Biocontainment certification by Public Health England, Porton Down, UK.

¹Maximum speed based on top performance models. Labware not shown to scale.

*Refers to the Legacy model.

Swinging buckets rotors

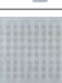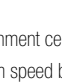

Thermo Scientific Rotor Name		BIOShield 720	TX-200	HIGHPlate 6000	Fiberlite H3-LV	M-20
Maximum Capacity (mL)		4 x 180 mL	4 x 180 mL	10 Microplates	28 Microplates	6 Microplates
Maximum Speed (rpm ¹)		6300 rpm	5500 rpm	6300 rpm	3600 rpm	4000 rpm
Maximum RCF (x g ¹)		7188 x g	5580 x g	6168 x g	2738 x g	2272 x g
Biocontainment 		✓ (Integrated)	✓ (Optional)	✓ (Integrated)		✓ (Optional)
Multifuge X4 Pro/Legend XF*				✓	✓	✓
Megafuge ST4 Plus						✓
Multifuge X1		✓	✓			✓
Megafuge ST16		✓	✓			✓
Maximum Capacity		1000 mL Bio-Bottle	-	-	-	-
		750 mL Bio-Bottle	-	-	-	-
		500 mL Nalgene Bottle	-	-	-	-
		500 mL Corning Bottle	-	-	-	-
		250 mL Nalgene Bottle	-	-	-	-
		250 mL Corning Bottle	-	-	-	-
		225 mL - 175 mL Conical Bottle	-	-	-	-
		50 mL Conical Tube	8	4	-	-
		50 mL Skirted Conical Tube	-	4	-	-
		50 mL Nalgene Oak Ridge Tube	-	-	-	-
		30 mL Sterilin Universal Tube	16	4	-	-
		15 mL Conical Tube	24	20	-	-
		10 mL Blood Collection Tube (16 x 100 mm)	48	28	-	-
		5/7 mL Blood Collection Tube (13 x 75-100 mm)	48	32	-	-
		4.5/6 mL Blood Collection Tube	-	32	-	-
		1.5/2 mL Microtube	160	48	-	-
		Microplate	-	10	28	6



Biocontainment certification by Public Health England, Porton Down, UK.

¹Maximum speed based on top performance models.
Labware not shown to scale.

*Refers to the Legacy model.

Rotors and adapters for swinging buckets rotors




TX-1000 Swinging Buckets Rotor¹, 90°, R_{max} 209 mm

Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4R Pro/Legend XFR*	4200	4122
Multifuge X4 Pro/Legend XF*	3800	3374
Megafuge ST4R Plus	4200	4122
Megafuge ST4 Plus	3800	3374
Multifuge X1	n/a	n/a
Megafuge ST16	n/a	n/a

¹Rotor, buckets and lids sold separately

Adapters

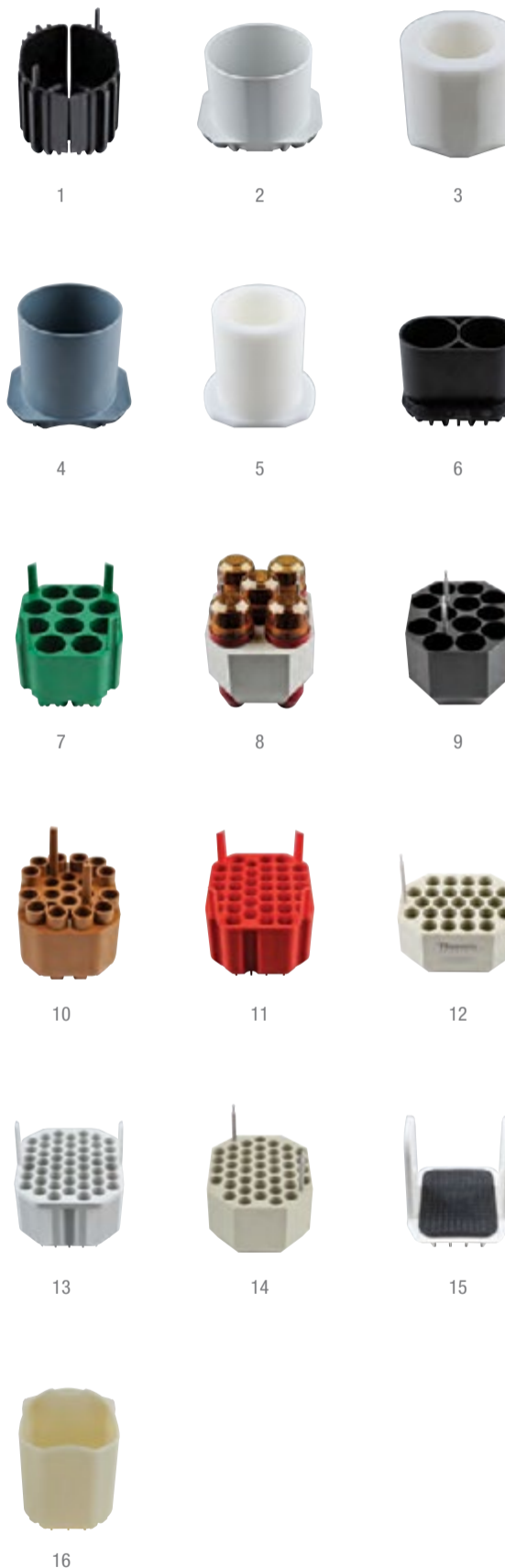

Cat. No.	Adapter Image	Description	Rotor Capacity (places x volume, mL)	Maximum Tube Dimensions (Ø x L, mm)
Rotor and Accessories				
75003017		TX-1000 Rotor Cross	4 x 1000 mL or 24 microplate	
75003001		TX-1000 Buckets (Set of 4)		
75007309		TX-1000 ClickSeal Biocontainment Lids (Set of 4)		
75007001		Replacement O-rings		
75007300		1000 mL Bio-Bottle - Polypropylene (Set of 4)		
Adapters for TX-1000 Rotor (Sets of 4)				
75007301	1	1000 mL Bio-Bottle	4 x 1000 mL	126 x 140
75007304	2	750 mL Polypropylene Bio-Bottle and any TX-750 Round adapter	4 x 750 mL	98 x 133
75004253	3	500 mL Nalgene Bottle	4 x 500 mL	70 x 160
75007302	4	500 mL Corning Bottle (requires Corning #431124)	4 x 500 mL	98 x 133
75005392	5	250 mL Corning Conical Bottle or 200 mL Nunc™ Bottle or 175 mL Nalgene Conical Bottle	4 x 250 mL	62 x 145
75007305	6	250 mL Nalgene Bottle	8 x 250 mL	62 x 145
75007305	6	225 mL BD Falcon™ (requires BD #352090)	8 x 225 mL	62 x 130
75007305	6	200 mL Nunc Conical Bottle (requires Nunc #377585)	8 x 200 mL	62.5 x 139
75007305	6	175 mL Nalgene Conical Bottle (requires Nalgene #DS3126-0175)	8 x 175 mL	62 x 130
75003674	7	50 mL Conical Tube	40 x 50 mL	30 x 120
75004255	8	50 mL Double Biocontainment Vessel for 50 mL Conical Tube (can be combined with ClickSeal lids)	20 x 50 mL	30 x 120
75004252	9	50 mL Nalgene Oak Ridge Tube	48 x 50 mL	29.5 x 120
75007306	10	15 mL Conical Tube	96 x 15 mL	17.5 x 121
75003672	11	10 mL Blood Collection Tube (16 x 100 mm) or Corex™/Kimble™ Tubes	148 x 10 mL	17 x 110
75003697	12	9/10 mL Blood Collection Tube (Sarstedt™)	100 x 9/10 mL	17 x 110
75003671	13	5/7 mL Blood Collection Tube (13 x 75-100 mm)	196 x 5/7 mL	14 x 110
75003709	14	4.5/6 mL Blood Collection Tube (Greiner™)	164 x 4.5/6 mL	14 x 110
75007303	15	Microplate Carriers	4 x 6 Standard or 4 x 2 Deep-Well	Height < 88 mm
75003829	16	Small Bag/Cell Culture Bags 4 x 2 Bags (< 350 mL)	4 x 2 Bags	

Biocontainment certification by Public Health England, Porton Down, UK.

*Refers to the Legacy model.

Rotors and adapters for swinging buckets rotors



TX-750 Swinging Buckets Rotor¹, 90°, R_{max} 195 mm

Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4 Pro/Legend XF*	4700	4816
Megafuge ST4 Plus	4700	4816
Multifuge X1	n/a	n/a
Megafuge ST16	n/a	n/a

¹Rotor, buckets and lids sold separately

Adapters



Cat. No.	Adapter Image	Description	Rotor Capacity (places x volume, mL)	Maximum Tube Dimensions (Ø x L, mm)
Rotor and Accessories				
75003180		TX-750 Rotor Cross		
75003608		TX-750 Round Buckets (Set of 4)	4 x 750 mL	
75003609		TX-750 Round ClickSeal Biocontainment Lids (Set of 4)		
75003610		Replacement TX-750 Round O-rings for lids (Set of 4)		
75006443		750 mL Bio-Bottle - Polypropylene (1 each)		
75003795	1	Microplate Carriers and T-75 Flask (includes tray and pads) (Set of 2)	4 Standard Microplate or 1 Deep well plate per carrier	Maximum height 60 mm
75003617	1	Microplate Carriers and T-75 Flask (includes tray and pads) (Set of 4)		
Adapters for TX-750 Rotor (Sets of 4, unless stated otherwise)				
Direct		750 mL Bio-Bottle	4 x 750 mL	98 x 133
Direct		500 mL Corning Conical Bottle (requires Corning #431124)	4 x 500 mL	98 x 133
75003792	2	250 mL Corning Conical Bottle (unsealed buckets only)	4 x 250 mL	62.5 x 139
75003710	3	250 mL Nalgene Bottle	4 x 250 mL	62.5 x 139
75003710	3	225 mL/175 mL BD Falcon Conical Bottle (requires BD #352090) (Open buckets only)	4 x 225 mL	62.5 x 139
75003710	3	200 mL Nunc Conical Bottle (requires Nunc #377585) (Open buckets only)	4 x 200 mL	62.5 x 139
75003710	3	175 mL Nalgene Conical Bottle (requires Nalgene #DS3126-0175) (Open buckets only)	4 x 175 mL	62.5 x 139
75003713	4	100 mL Round Bottom Open-Top Tube	8 x 100 mL	45 x 123
75003714	5	50 mL Conical Tube (includes sealing vessel) can be combined with ClickSeal lids	12 x 50 mL	30 x 120
75003638	6	50 mL Conical Tube	28 x 50 mL	30 x 120
75003824	7	50 mL Conical or Skirted Tube	20 x 50 mL	30 x 120
75003715	8	50 mL Nalgene Oak Ridge Tube	28 x 50 mL	29.5 x 120
75003716	9	30 mL Sterilin Universal Tube	28 x 30 mL	25.5 x 120
75003639	10	15 mL Conical Tube	56 x 15 mL	17 x 122
75003719	11	15 mL Blood Collection Tube (17 x 125 mm) (inner circle only)	32 x 15 mL	17 x 134
75003719	11	10 mL Blood Collection (16 x 100 mm) or 15 mL Corex/Kimble Tube	84 x 10/15 mL	17 x 117
75003718	12	14 mL Round or Conical Urine Tube	48 x 14 mL	18 x 124
75003723	13	5/7 mL or 4.5/6 mL Blood Collection Tube (13 x 75-100 mm)	108 x 5/7 mL	13 x 116
75003724	14	5 mL RIA or Round Bottom Tube (without cap)	148 x 5 mL	13 x 111
75003732	14	5/7 mL Round Bottom Tube (without cap) with Decanting Aid	148 x 5/7 mL	13 x 111
75003733	15	1.5/2 mL Microtube	192 x 2 mL	11 x 50
75008383	16	T-75 Nunc Easy Flask (Set of 2)	4 x T-75	-
75008384	17	T-25 Nunc Easy Flask (Set of 2)	8 x T-25	-

Biocontainment certification by Public Health England, Porton Down, UK.

*Refers to the Legacy model.

Rotors and adapters for swinging buckets rotors



TX-750 Swinging Buckets Rotor, 90°, R_{max} 195 mm

Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4R Pro/Legend XFR*, Megafuge ST4R Plus, 230V	4700	4816
Multifuge X4 Pro/Legend XF*, Megafuge ST4 Plus, 230V	4500	4415
Multifuge X4 Pro Series/Legend XF*, Megafuge ST4 Plus Series, 120V	4300	4031
Multifuge X1 or Megafuge ST16	n/a	n/a

Adapters



1 2 3

Cat. No.	Adapter Image	Description	Rotor Capacity (places x volume, mL)	Maximum Tube Dimensions (Ø x L, mm)
Rotor and Accessories				
75003180		TX-750 Rotor Cross		
75003614		TX-750 Rectangular Buckets (Set of 4)	4 x 250 mL	
75003615		TX-750 Rectangular Click Seal Biocontainment Lids (Set of 4)		
75003616		Replacement TX-750 Rectangular O-rings for lids (Set of 4)		
Adapters for TX-750 Rotor (Sets of 4)				
75003737	1	250 mL Flat Bottom Bottle	4 x 250 mL	63 x 130
75003738	2	150 mL Round Bottom Open-Top Tube	4 x 150 mL	56 x 124
75003742	3	100 mL Round Bottom Open-Top Tube	4 x 100 mL	46 x 124
75003685	4	50 mL Conical Tube	20 x 50 mL	30 x 121
75003749	5	50 mL Round Bottom Tube	12 x 50 mL	35 x 121
75003750	6	45 mL Flat/Round Tube	16 x 45 mL	30 x 121
75003755	7	30 mL Sterilin Universal Tube	16 x 30 mL	26 x 121
75003756	8	25 mL DIN Round Bottom Tube	24 x 25 mL	26 x 121
75003684	9	15 mL Conical Tube	48 x 15 mL	17 x 120
75003758	10	14 mL Flanged Round Bottom Tube	36 x 14 mL	18 x 121
75003759	11	14 mL Round or Conical Urine Tube	36 x 15 mL	17 x 124
75003767	12	10 mL Blood Collection (16 x 100 mm) or 15 mL DIN Tube	64 x 10/15 mL	17 x 118
75003768	13	5/7 mL or 4.5/6 mL Blood Collection Tube	80 x 5/7 mL	14 x 118
75003769	14	5/7 mL Round Bottom Tube (13 x 75-100 mm)	112 x 5/7 mL	13 x 119
75003770	15	1.5/2 mL Microtube	224 x 2 mL	11 x 50



4 5 6



7 8 9



10 11 12



BIOShield 1000A Swinging Buckets Rotor, 90°, R_{max} 178 mm

Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4 Pro/Legend XF*	6000	7164
Megafuge ST4 Plus	5300	5590
Multifuge X1	n/a	n/a
Megafuge ST16	n/a	n/a

Cat. No.	Adapter Image	Description	Rotor Capacity (places x volume, mL)	Maximum Tube Dimensions (Ø x L, mm)
Rotor and Accessories				
75003182		BIOShield 1000A	4 x 250 mL	
Adapters for BIOShield 1000A Rotor (Sets of 4)				
75003737	1	250 mL Flat Bottom Bottle	4 x 250 mL	63 x 130
75003738	2	150 mL Round Bottom Open-Top Tube	4 x 150 mL	56 x 124
75003742	3	100 mL Round Bottom Open-Top Tube	4 x 100 mL	46 x 124
75003749	5	50 mL Round Bottom Tube	12 x 50 mL	35 x 121
75003750	6	45 mL Flat/Round Tube	16 x 45 mL	30 x 121
75003755	7	30 mL Sterilin Universal Tube	16 x 30 mL	30 x 121
75003756	8	25 mL DIN Round Bottom Tube	24 x 25 mL	26 x 121
75003758	10	14 mL Flanged Round Bottom Tube	36 x 14 mL	18 x 121
75003759	11	14 mL Round or Conical Urine Tube	36 x 15 mL	17 x 124
75003767	12	10 mL Blood Collection (16 x 100 mm) or 15 mL DIN Tube	64 x 10/15 mL	17 x 118
75003768	13	5/7 mL or 4.5/6 mL Blood Collection Tube	80 x 5/7 mL	14 x 118
75003769	14	5/7 mL Round Bottom Tube (13 x 75-100 mm)	112 x 5/7 mL	13 x 119
75003770	15	1.5/2 mL Microtube	224 x 2 mL	11 x 50
75003643	16	50 mL Conical Tube	12 x 50 mL	30 x 121
75003642	17	15 mL Conical Tube	36 x 15 mL	17 x 120



13 14 15



16 17

Biocontainment certification by Public Health England, Porton Down, UK.

*Refers to the Legacy model.

Rotors and adapters for swinging buckets rotors




TX-400 Swinging Buckets Rotor¹, 90°, R_{max} 168 mm

Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4 Pro/Legend XF*	n/a	n/a
Megafuge ST4 Plus	n/a	n/a
Multifuge X1	5000	4696
Megafuge ST16	5000	4696

¹Rotor, buckets and lids sold separately

Adapters



Cat. No.	Adapter Image	Description	Rotor Capacity (places x volume, mL)	Maximum Tube Dimensions (Ø x L, mm)
Rotor and Accessories				
75003181		TX-400 Rotor Cross		
75003655		TX-400 Round Buckets (Set of 4)	4 x 400 mL	
75003656		TX-400 Round ClickSeal Biocontainment Lids (Set of 4) 		
75003657		Replacement TX-400 O-rings for lids (Set of 4)		
75007585		400 mL Bio-Bottle - Polypropylene (Set of 12)		
Adapters for TX-400 Rotor (Sets of 4)				
Direct		400 mL Bio-Bottle	4 x 400 mL	80 x 125
75003788	1	250 mL Thermo Scientific Nalgene™ Bottle	4 x 250 mL	62 x 135
75003788	1	200 mL Thermo Scientific Nunc Conical Bottle (requires Nunc #377585)	4 x 200 mL	62 x 125
75003788	1	225 mL/175 mL BD Falcon Conical Bottle (requires BD #352090)	4 x 175 mL	62 x 125
75003708	2	100 mL Round Bottom Open-Top Tube	4 x 100 mL	45 x 117
75003707	3	50 mL DIN Round Bottom Tube	12 x 50 mL	34.5 x 105
75003799	4	50 mL Nalgene™ Oak Ridge Tube	16 x 50 mL	29.5 x 114
75003683	5	50 mL Conical Tube	16 x 50 mL	30 x 116
75003703	6	30/25 mL DIN Round/Flat Bottom Tube	20 x 30 mL	25.5 x 108
75003706	7	30 mL Sterilin Universal Tube	12 x 30 mL	25 x 110
75003682	8	15 mL Conical Tube	36 x 15 mL	17 x 121
75003794	9	15 mL Blood Collection Tube (17 x 125 mm)	16 x 15 mL	17 x 131
75003704	10	15 mL Round Bottom Tube (Sarstedt)	40 x 15 mL	17 x 105
75003798	11	14 mL Round or Conical Urine Tube	28 x 14 mL	17 x 112
75003681	12	10 mL Blood Collection (16 x 100 mm)	56 x 10/15 mL	16 x 113
75003680	13	5/7 mL Blood Collection Tube (13 x 75-100 mm)	76 x 5/7 mL	13 x 110
75003825	14	4.5/6 mL Blood Collection Tube (Greiner)	64 x 4.5/6 mL	13 x 110
75003793	15	3 mL RIA or Round Bottom Tube (without cap)	76 x 3/5 mL	11 x 110
75003700	16	1.5/2 mL Microtube	136 x 2 mL	11 x 45

 Biocontainment certification by Public Health England, Porton Down, UK.

*Refers to the Legacy model.


Rotors and adapters for swinging buckets rotors



BIOShield 720 Swinging Buckets Rotor, 90°, R_{max} 162 mm

Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4 Pro/Legend XF*	n/a	n/a
Megafuge ST4 Plus	n/a	n/a
Multifuge X1	6300	7188
Megafuge ST16	5300	5088

Adapters

Cat. No.	Adapter Image	Description	Rotor Capacity (places x volume, mL)	Maximum Tube Dimensions (Ø x L, mm)
Rotor and Accessories				
75003183		BIOShield 720	4 x 180 mL	n/a
75003693		180 mL Bio-Bottle - Polypropylene (Set of 12)		
BIOShield 720A Adapters (Sets of 4)				
75003813	1	150 mL Round Bottom Open-Top Tube	4 x 150 mL	50 x 122
75003814	2	100 mL Round Bottom Open-Top Tube	4 x 100 mL	44 x 120
75003677	3	50 mL Conical Tube	8 x 50 mL	30 x 116
75003816	4	50 mL DIN Round Bottom Tube	8 x 50 mL	34 x 109
75003818	5	30 mL Sterilin Universal Tube	16 x 30 mL	26 x 105
75003817	6	25 mL DIN Round Bottom Tube	20 x 25 mL	25 x 107
75003678	7	15 mL Conical Tube	24 x 15 mL	17 x 120
75003820	8	15 mL Sarstedt Tube	44 x 15 mL	17 x 103
75003701	9	10 mL Blood Collection Tube (16 x 100 mm)	48 x 10 mL	17 x 105
75003821	10	5/7 mL Blood Collection Tube (13 x 75-100 mm)	48 x 5/7 mL	13 x 106
75003822	11	5/7 mL Round Bottom Open-Top Tube	80 x 5/7 mL	13 x 104
75003823	12	1.5/2 mL Conical Microtube	160 x 2 mL	11 x 45



1



2



3



4



5



6



7



8



9



10



11



12



Biocontainment certification by Public Health England, Porton Down, UK.

*Refers to the Legacy model.

Rotors and adapters for swinging buckets rotors




TX-200 Swinging Buckets Rotor¹, 90°, R_{max} 165 mm

Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4 Pro/Legend XF*	n/a	n/a
Megafuge ST4 Plus	n/a	n/a
Multifuge X1R 230V and 120V	5500	5580
Multifuge X1 120V	5250	5084
Megafuge ST16 R 230V and 120V	5500	5580
Megafuge ST16 120V	5250	5084

¹Rotor, buckets and lids sold separately

Adapters



Cat. No.	Adapter Image	Description	Rotor Capacity (places x volume, mL)	Maximum Tube Dimensions (Ø x L, mm)
Rotor and Accessories				
75003658		TX-200 Rotor Cross		
75003659		TX-200 Round Buckets (Set of 4)	4 x 180 mL	
75003660		TX-200 Round ClickSeal Biocontainment Lids (Set of 4) 		
75003687		Replacement TX-200 O-rings for lids (Set of 4)		
75003800		180 mL Bio-Bottle - Polypropylene (Set of 12)		
Adapters for TX-200 Rotor (Sets of 4)				
Direct		180 mL Bio-Bottle	4 x 180 mL	56.5 x 123
75003801	1	100 mL Round Bottom Open-Top Tube	4 x 100 mL	48 x 122
75003802	2	50 mL DIN Round Bottom Tube	4 x 50 mL	35 x 126
75003815	3	50 mL Round Bottom Tube	4 x 50 mL	32 x 126
75003803	4	50 mL Conical or Skirted Tube	4 x 50 mL	30 x 126
75003804	5	30 mL Sterilin™ Universal Tube	4 x 30 mL	26 x 118
75003805	6	25 mL DIN Round/Flat Bottom Tube	8 x 25 mL	25 x 118
75003806	7	20 mL Round Bottom Tube	12 x 25 mL	22 x 125
75003771	8	15 mL Conical Tube	20 x 20 mL	17 x 121
75003809	9	15 mL Blood Collection Tube	16 x 15 mL	17 x 121
75003807	10	14 mL Conical Urine Tube	8 x 14 mL	17 x 114
75003808	11	10 mL Blood Collection (16 x 100 mm) or 15 mL Corex/Kimble Tube	28 x 10/15 mL	17 x 121
75003785	12	5/7 mL or 4.5/6 mL Blood Collection Tube (13 x 75-100 mm)	32 x 5/7 mL	13 x 122
75003810	13	5/7 mL Round Bottom Open-Top Tube	36 x 5/7 mL	13 x 118
75003811	14	3/5 mL RIA or Round Bottom Tube (without cap)	48 x 3/5 mL	12 x 124
75003812	15	1.5/2 mL Microtube	48 x 2 mL	11 x 45

 Biocontainment certification by Public Health England, Porton Down, UK.


*Refers to the Legacy model.

Rotors and adapters for swinging buckets rotors



HIGHPlate 6000 Swinging Buckets Rotor, 90°, R_{max} 139 mm

Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4 Pro/Legend XF*	6300	6168
Megafuge ST4 Plus	n/a	n/a
Multifuge X1	n/a	n/a
Megafuge ST16	n/a	n/a

Cat. No.	Description	Rotor Capacity (places x volume, mL)	Maximum Height
Rotor			
75003606	HIGHPlate 6000 Rotor 	2 x 5 Standard Microplate or 2 x 1 Filter plates	< 80 mm



Fiberlite H3-LV Swinging Buckets Rotor, 90°, R_{max} 189 mm

Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4 Pro/Legend XF*	3600	2738
Megafuge ST4 Plus	n/a	n/a
Multifuge X1	n/a	n/a
Megafuge ST16	n/a	n/a


Cat. No.	Description	Rotor Capacity (places x volume, mL)	Maximum Height
Rotor			
75003665	Fiberlite H3-LV Rotor	4 x 7 Standard Microplate or 2 x 4 Deepwell	< 100 mm



M-20 Swinging Buckets Rotor¹, 90°, R_{max} 127 mm

Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4 Pro/Legend XF*	4000	2272
Megafuge ST4 Plus	4000	2272
Multifuge X1	4000	2272
Megafuge ST16	4000	2272

¹ Biocontainment system sold separately



Cat. No.	Description	Rotor Capacity (places x volume, mL)	Maximum Height
Rotor and Accessories			
75003624	M-20 Microplate Rotor including standard open carriers	2 x 3 Standard Microplate 2 x 1 Deepwell	< 65 mm
75003625	ClickSeal Carriers including Carrier and lids 		< 45 mm
75002011	ClickSeal Replacement lids (Set of 2)		
75002012	Replacement M-20 O-rings for lids		

 Biocontainment certification by Public Health England, Porton Down, UK.

*Refers to the Legacy model.

Fixed angle rotors



Thermo Scientific Rotor Name		Fiberlite F14-6x250 LE	Fiberlite F15-6x100y	Fiberlite F10-6x100 LEX	Fiberlite F13-14x50cy	Fiberlite F15-8x50cy	CLINIConic
Maximum Capacity (mL)		6 x 250 mL	6 x 100 mL	6 x 100 mL	14 x 50 mL	8 x 50 mL	30 x 15 mL
Maximum Speed (rpm ¹)		11000 rpm	15000 rpm	10000 rpm	10500 rpm	14500 rpm	5650 rpm
Maximum RCF (x g ¹)		18533 x g	24652 x g	17105 x g	15038 x g	24446 x g	4997 x g
Biocontainment 		✓	✓	✓	✓	✓	
Multifuge X4 Pro/Legend XF*		✓	✓	✓	✓	✓	
Megafuge ST4 Plus			✓				
Multifuge X1			✓	✓	✓	✓	✓
Megafuge ST16			✓				✓
Maximum Capacity	 250 mL Nalgene Bottle	6	–	–	–	–	–
	 85 mL Nalgene Oak Ridge Tube	6	6	6	–	–	–
	 50 mL Nalgene Oak Ridge Tube	6	6	6	14	8	–
	 50 mL Conical Tube	6	6	6	14	8	–
	 30 mL Nalgene Oak Ridge Tube	12	6	6	14	8	–
	 16 mL Nalgene Oak Ridge Tube	30	6	6	14	8	–
	 15 mL Conical Tube	30	6	6	14	8	30
	 10 mL Nalgene Oak Ridge Tube	42	12	12	14	8	–
	 10 mL Blood Collection Tube (16 x 100 mm)	–	–	–	–	–	30
	 5/7 mL Blood Collection Tube (13 x 75-100 mm)	–	–	–	–	–	30
	 6.5 mL Round Tube	–	12	12	–	–	–
	 1.5/2 mL Microtube	–	24	24	–	–	–
















Biocontainment certification by Public Health England, Porton Down, UK.

¹ Maximum speed based on top performance models. Labware not shown to scale.

*Refers to the Legacy model.

Fixed angle rotors



Thermo Scientific Rotor Name	8 x 50 Sealed	HIGHConic II	MicroClick 18 x 5	Fiberlite F21-48x2	Microliter 48 x 2	MicroClick 30 x 2	Microliter 30 x 2
Maximum Capacity (mL)	8 x 50 mL	6 x 100 mL	18 x 5 mL	48 x 2 mL	48 x 2 mL	30 x 2 mL	30 x 2 mL
Maximum Speed (rpm ¹)	6700 rpm	10350 rpm	15000 rpm	15200 rpm	15200 rpm	14000 rpm	15200 rpm
Maximum RCF (x g ¹)	7177 x g	15090 x g	24652 x g	25055 x g	25314 x g	21693 x g	25830 x g
Biocontainment 	✓ (Individual)	✓	✓	✓	✓	✓	✓
Multifuge X4 Pro/Legend XF*		✓	✓	✓	✓	✓	✓
Megafuge ST4 Plus		✓	✓	✓	✓	✓	✓
Multifuge X1	✓	✓	✓	✓	✓	✓	✓
Megafuge ST16	✓	✓	✓	✓	✓	✓	✓
Maximum Capacity 	250 mL Nalgene Bottle	-	-	-	-	-	-
	85 mL Nalgene Oak Ridge Tube	-	6	-	-	-	-
	50 mL Nalgene Oak Ridge Tube	-	6	-	-	-	-
	50 mL Conical Tube	8	6	-	-	-	-
	30 mL Nalgene Oak Ridge Tube	-	6	-	-	-	-
	16 mL Nalgene Oak Ridge Tube	-	6	-	-	-	-
	15 mL Conical Tube	8	6	-	-	-	-
	10 mL Nalgene Oak Ridge Tube	-	12	-	-	-	-
	10 mL Blood Collection Tube (16 x 100 mm)	8	-	-	-	-	-
	5/7 mL Blood Collection Tube (13 x 75-100 mm)	8	-	-	-	-	-
	6.5 mL Round Tube	-	12	-	-	-	-
	5 mL microtube/conical	-	-	18	-	-	-
	1.5/2 mL Microtube	-	24	18	48	48	30 (Including spin columns)

 Biocontainment certification by Public Health England, Porton Down, UK.

¹Maximum speed based on top performance models. Labware not shown to scale.

*Refers to the Legacy model.

Rotors and adapters for fixed angle rotors



Fiberlite F14-6x250 LE Fixed Angle Rotor, 23°, R_{max} 137 mm

Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4 Pro/Legend XF* 230V	11000	18533
Multifuge X4 Pro/Legend XF* 120V	10000	15317
Megafuge ST4 Plus	n/a	n/a
Multifuge X1	n/a	n/a
Megafuge ST16	n/a	n/a

Cat. No.	Description	Rotor Capacity (places x volume, mL)	Maximum Tube Dimensions (Ø x L, mm)
Rotor and Accessories			
75003662	Fiberlite F14-6 x 250LE Rotor for 250 mL Nalgene Bottles	6 x 250 mL	62 x 125
021-062153	Replacement O-ring Kit		
3141-0250	250 mL Nalgene Polypropylene Bottle (Set of 36)		
3140-0251	250 mL Nalgene Polycarbonate Bottle (Set of 36)		
Adapters for F14-6 X 250LE Rotor (Sets of 2, unless stated otherwise)			
010-1119	85 mL Nalgene Oak Ridge Tube	6 x 85 mL	38 x 109
010-0138	50 mL Nalgene Oak Ridge Tube	6 x 50 mL	29.5 x 121
75100136	50 mL Conical Tube	6 x 50 mL	30 x 121
010-1072	30 mL Nalgene Oak Ridge Tube / 38 mL Round Tube	12 x 30 mL	26 x 105
010-1074	16 mL Nalgene Oak Ridge Tube	30 x 16 mL	18 x 107
010-1410	15 mL Conical Tube	30 x 15 mL	17 x 121
00389	10 mL Nalgene Oak Ridge Tube / 12 mL Round Tube (each)	42 x 10 mL	17 x 110



Fiberlite F10-6x100 LEX Fixed Angle Rotor, 45°, R_{max} 122 mm Fiberlite F15-6x100y Fixed Angle Rotor, 25°, R_{max} 98 mm

Centrifuge	F10 Maximum Speed (rpm)	F10 Maximum RCF (x g)	F15 Maximum Speed (rpm)	F15 Maximum RCF (x g)
Multifuge X4 Pro/Legend XF*	10500	15038	15000	24652
Megafuge ST4 Plus	10500	15038	13000	18516
Multifuge X1	10500	15038	15000	24652
Megafuge ST16	10500	15038	13000	18516

Cat. No.	Description	Rotor Capacity (places x volume, mL)	Maximum Tube Dimensions (Ø x L, mm)
Rotor and Accessories			
75003698	Fiberlite F15-6 x 100y Rotor for 100 mL/85 mL Nalgene	6 x 100 mL	38 x 106
75003340	Fiberlite F10-6 x 100 LEX Rotor for 100mL/85mL Nalgene	6 x 100 mL	38 x 106
021-069031	Replacement O-ring Kit		
3118-0085	85 mL Nalgene Oak Ridge Tube (Set of 50)	6 x 85 mL	38 x 106
Adapters for F15-6 x 100y Rotor (Each, unless stated otherwise)			
75003102	50 mL Nalgene Oak Ridge Tube	6 x 50 mL	29.5 x 117
75003103	50 mL Conical Tube	6 x 50 mL	30 x 121
75003094	30 mL Nalgene™ Oak Ridge Tube / 38 mL Round Tube	6 x 30 mL	26 x 110
76002906	16 mL Nalgene Oak Ridge Tube	6 x 16 mL	18 x 123
75003095	15 mL Conical Tube	6 x 15 mL	17 x 123
75003093	10 mL Nalgene Oak Ridge Tube / 12 mL Round Tube	12 x 12 mL	17 x 95
75003092	6.5 mL Round Bottom Tube	12 x 6.5 mL	13 x 114
75003091	1.5/2 mL Microtube	24 x 2 mL	11 x 40



Biocontainment certification by Public Health England, Porton Down, UK.


*Refers to the Legacy model.

Rotors and adapters for fixed angle rotors



Fiberlite F13-14x50cy Fixed Angle Rotor, 34°, R_{max} 153 mm


Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4 Pro/Legend XF* 230V	10000	17105
Multifuge X4 Pro/Legend XF* 120V	9250	14636
Megafuge ST4 Plus	n/a	n/a
Multifuge X1	8500	12359
Megafuge ST16	n/a	n/a

Cat. No.	Description	Rotor Capacity (places x volume, mL)	Maximum Tube Dimensions (Ø x L, mm)
Rotor and Accessories			
75003661	Fiberlite F13-14x50cy for 50 mL Conical Tubes 	14 x 50 mL	30 x 121
021-149027	Replacement O-ring Kit		
Adapters for F13-14x50cy Rotor (Set of 2, unless otherwise stated)			
010-0377	50 mL Nalgene Oak Ridge Tube	14 x 50 mL	29.5 x 108
010-1147	30 mL Nalgene Oak Ridge Tube	14 x 30 mL	26 x 102
010-0376	16 mL Nalgene Oak Ridge Tube	14 x 16 mL	18 x 123
75100378	15 mL Conical Tube	14 x 15 mL	17 x 121
010-1340	15 mL Nunc Conical Tube and Filtration Devices	14 x 15 mL	17 x 121
010-1311	10 mL Nalgene Oak Ridge Tube / 12 mL Round Tube	28 x 12 mL	17 x 95



Fiberlite F15-8x50cy Fixed Angle Rotor, 25°, R_{max} 104 mm

Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4 Pro/Legend XF*	14500	24446
Megafuge ST4 Plus	n/a	n/a
Multifuge X1 230V	14500	24446
Multifuge X1 120V	14000	22789
Megafuge ST16	n/a	n/a

Cat. No.	Description	Rotor Capacity (places x volume, mL)	Maximum Tube Dimensions (Ø x L, mm)
Rotor and Accessories			
75003663	Fiberlite F15-8 x 50cy for 50 mL Conical Tubes 	8 x 50 mL	30 x 121
021-085077	Replacement O-ring Kit		
Adapters for F15-8x50cy Rotor (Set of 2, unless otherwise stated)			
010-0377	50 mL Nalgene Oak Ridge Tube	8 x 50 mL	29 x 108
010-1147	30 mL Nalgene Oak Ridge Tube	8 x 30 mL	26 x 102
010-0376	16 mL Nalgene Oak Ridge Tube	8 x 16 mL	18 x 123
75100378	15 mL Conical Tube	8 x 15 mL	17 x 121
010-1340	15 mL Nunc Conical Tube and Filtration Devices	8 x 15 mL	17 x 121
010-1311	10 mL Nalgene Oak Ridge Tube / 12 mL Round Tube	8 x 12 mL	17 x 95



Biocontainment certification by Public Health England, Porton Down, UK.

*Refers to the Legacy model.

Rotors and adapters for fixed angle rotors



CLINIConic Fixed Angle Rotor, 37°, R_{max} 140 mm


Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4 Pro/Legend XF*	n/a	n/a
Megafuge ST4 Plus	n/a	n/a
Multifuge X1	5650	4997
Megafuge ST16	5650	4997

Cat. No.	Description	Rotor Capacity (places x volume, mL)	Maximum Tube Dimensions (Ø x L, mm)
Rotor			
75003623	CLINIConic Fixed Angle Rotor for 15 mL Round and Conical Tubes	30 x 15 mL	17 x 121
Adapters for CLINIConic Rotor (Each)			
75003702	10 mL Round Bottom Tube	30 x 10 mL	17 x 131
11172596	7 mL Blood Collection Tubes (13 x 100 mm)	30 x 7 mL	17 x 131
11172595	5 mL Blood Collection Tubes (13 x 75 mm)	30 x 5 mL	13 x 95



8 x 50 Individually Sealed Fixed Angle Rotor, 45°, R_{max} 143 mm

Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4 Pro/Legend XF*	n/a	n/a
Megafuge ST4 Plus	n/a	n/a
Multifuge X1 230V	6700	7177
Multifuge X1 120V	6500	6755
Megafuge ST16 230V	6700	7177
Megafuge ST16 120V	6500	6755

Cat. No.	Description	Rotor Capacity (places x volume, mL)	Maximum Tube Dimensions (Ø x L, mm)
Rotor and Accessories			
75003694	8 x 50 Individually Sealed for 50 mL Conical Tubes 	8 x 50 mL	30 x 121
75003689	Replacement O-ring Kit		
Adapters for 8 x 50 Sealed Rotor (Set of 6)			
75005755	15 mL Conical Tube	8 x 50 mL	17 x 121
75005747	10 mL Blood Collection Tube (16 x 100 mm)	8 x 10 mL	16 x 110
75005748	7 mL Blood Collection Tube (13 x 100 mm)	8 x 7 mL	13 x 110
75005749	3.5 mL Blood Collection Tube	16 x 3.5 mL	11 x 75



Biocontainment certification by Public Health England, Porton Down, UK.


*Refers to the Legacy model.

Rotors and adapters for fixed angle rotors



HIGHConic II Aluminium Fixed Angle Rotor, 45°, R_{max} 126 mm


Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4 Pro/Legend XF*	10350	15090
Megafuge ST4 Plus	10350	15090
Multifuge X1	10350	15090
Megafuge ST16	10350	15090

Cat. No.	Description	Rotor Capacity (places x volume, mL)	Maximum Tube Dimensions (Ø x L, mm)
Rotor and Accessories			
75003620	HIGHConic Fixed Angle Rotor for 100 mL/85 mL Nalgene Oak Ridge Tube 	6 x 100 mL	38 x 123
75003058	Replacement O-ring Kit		
3118-0085	85 mL Nalgene Oak Ridge Tube (Set of 50)	6 x 85 mL	38 x 106
Adapters for HIGHConic Rotor (Each, unless stated)			
75003102	50 mL Nalgene Oak Ridge Tube	6 x 50 mL	30 x 117
75003103	50 mL Conical Tube	6 x 50 mL	30 x 121
75003094	30 mL Nalgene Oak Ridge Tube / 38 mL Round Tube	6 x 30 mL	26 x 110
76002906	16 mL Nalgene Oak Ridge Tube	6 x 16 mL	18 x 123
75003095	15 mL Conical Tube	6 x 15 mL	17 x 123
75003093	10 mL Nalgene Oak Ridge Tube / 12 mL Round Tube	12 x 12 mL	17 x 95
75003092	6.5 mL Round Bottom Tube	12 x 6.5 mL	13 x 114
75003091	1.5/2 mL Microtube	24 x 2 mL	11 x 40



MicroClick 18 x 5 Aluminum Fixed Angle Rotor, 45°, R_{max} 102 mm


Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4 Pro/Legend XF*	15000	24652
Megafuge ST4 Plus	15000	24652
Multifuge X1	15000	24652
Megafuge ST16	15000	24652

Cat. No.	Description	Rotor Capacity (places x volume, mL)	Maximum Tube Dimensions (Ø x L, mm)
Rotor and Accessories			
75005765	18 x 5 MicroClick Rotor for 5mL Conical and Flip Top Tubes 	18 x 5 mL	17 x 6
20059119	ClickSeal	n/a	n/a
75005726	O-Ring Set	n/a	n/a
75003500	O-Ring Grease	n/a	n/a
Adapters for MicroClick Rotor			
75005756	1.5/2 mL Microtube	18 x 5	11 x 45



Fiberlite F21-48 x 2 Fixed Angle Rotor, 45°, R_{max} 97 mm

Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4 Pro/Legend XF*	15200	25055
Megafuge ST4 Plus	15200	25055
Multifuge X1	15200	25055
Megafuge ST16	15200	25055

Cat. No.	Description	Rotor Capacity (places x volume, mL)	Maximum Tube Dimensions (Ø x L, mm)
Rotor and Accessories			
75003664	Fiberlite F21-48 x 2 for 1.5/2 mL Microtubes 	48 x 2 mL	11 x 50
021-489021	Replacement O-ring Kit		



Biocertification by Public Health England, Porton Down, UK.

*Refers to the Legacy model.

Rotors and adapters for fixed angle rotors



Microliter 48 x 2 Aluminum Fixed Angle Rotor, 45°, R_{max} 98 mm

Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4 Pro/Legend XF*	15200	25314
Megafuge ST4 Plus	15200	25314
Multifuge X1	15200	25314
Megafuge ST16	15200	25314

Cat. No.	Description	Rotor Capacity (places x volume, mL)	Maximum Tube Dimensions (Ø x L, mm)
Rotor and Accessories			
75003602	48 x 2 Aluminum for 1.5/2 mL Microtubes	48 x 2 mL	11 x 50
75003349	Replacement O-ring Kit		



MicroClick 30 x 2 Aluminum Fixed Angle Rotor, 45°, R_{max} 99 mm

Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4 Pro/Legend XF*	14000	21693
Megafuge ST4 Plus	14000	21693
Multifuge X1	14000	21693
Megafuge ST16	14000	21693

Cat. No.	Description	Rotor Capacity (places x volume, mL)	Maximum Tube Dimensions (Ø x L, mm)
Rotor and Accessories			
75005719	30 x 2 MicroClick for 1.5/2 mL Microtubes and Spin Columns	30 x 2 mL	11 x 50
75005726	Replacement O-ring Kit		



Microliter 30 x 2 Aluminum Fixed Angle Rotor, 45°, R_{max} 100 mm

Centrifuge	Maximum Speed (rpm)	Maximum RCF (x g)
Multifuge X4 Pro/Legend XF*	15200	25830
Megafuge ST4 Plus	15200	25830
Multifuge X1	15200	25830
Megafuge ST16	15200	25830

Cat. No.	Description	Rotor Capacity (places x volume, mL)	Maximum Tube Dimensions (Ø x L, mm)
Rotor and Accessories			
75003652	30 x 2 Aluminum for 1.5/2 mL Microtubes	30 x 2 mL	11 x 50
75003349	Replacement O-ring Kit		



Biocontainment certification by Public Health England, Porton Down, UK.

Microliter Rotor Adapters

Cat. No.	Description	Maximum Tube Dimensions (Ø x L, mm)
Adapters for all Microliter Rotors (Set of 30)		
75005753	0.5 mL Microtube	8 x 44
75005754	0.25 mL Microtube	6 x 46
75005752	0.2 mL PCR Tube	5 x 20

000 «Диаэм»

Москва
ул. Магаданская, д. 7, к. 3 ■ тел./факс: (495) 745-0508 ■ sales@dia-m.ru

www.dia-m.ru

С.-Петербург
+7 (812) 372-6040
spb@dia-m.ru

Новосибирск
+7 (383) 328-0048
nsk@dia-m.ru

Воронеж
+7 (473) 232-4412
vrn@dia-m.ru

Йошкар-Ола
+7 (927) 880-3676
nba@dia-m.ru

Красноярск
+7 (923) 303-0152
krsk@dia-m.ru

Казань
+7 (843) 210-2080
kazan@dia-m.ru

Ростов-на-Дону
+7 (863) 303-5500
rnd@dia-m.ru

Екатеринбург
+7 (912) 658-7606
ekb@dia-m.ru

Кемерово
+7 (923) 158-6753
kemerovo@dia-m.ru

Армения
+7 (094) 01-0173
armenia@dia-m.ru

